

TICKER® Unlimited & Gentle Dose®

Patient Instructions


Remove the safety cap on TICKER® Unlimited by placing your thumb under the thumb-tab and pushing upward.

PRIMING THE TICKER UNIT

If you begin to turn the applicator and the cream does not come out, the dispenser may need to be primed. **Priming:** Slowly turn the applicator base to the right until cream begins to dispense. (It may require a few turns) Stop when the cream appears at the top. Use a tissue to wipe off the excess. Your device is now primed and ready to use!

Attach the Gentle Dose® applicator to the nozzle of TICKER® Unlimited. Turn TICKER® Unlimited dial to the right to get your dosage as specified on the prescription. DIAL as many CLICKS as instructed. (Note the amount being dispensed will appear on the dosing segments of the plunger. Each segment is arranged in 0.5mL intervals.

The plunger of the applicator is removable.

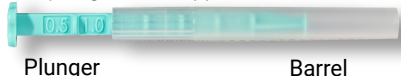

After dialing your dosage, wait at least 5-10 seconds in order to give the preparation ample time to exit through the nozzle and collect inside Gentle Dose®. The Gentle Dose® applicator is now ready to be inserted.

Follow your doctor or pharmacist's instructions for how to apply or insert the custom medication. Press the plunger all the way down to dispense the correct dose, then remove the applicator from the body. Sanitize the applicator in lukewarm running water and let it air dry.

When the top of the piston of the dispenser reaches the extruded word "REFILL," it is time to refill your medication. There are approximately 6-7 milliliters remaining and it may be time to contact your pharmacy for additional refills.

Thank You
for choosing our pharmacy

Gentle Dose® Cleaning Instructions

- Depress the plunger all the way to ensure all of the medication was dispensed
- Remove plunger and barrel from their counterparts by pulling them away from each other
- Wash plunger with lukewarm water. Use mild, hypoallergenic soap when desired. You may use a thin brush inside the barrel if necessary
- Allow the parts to air dry. You may then connect all parts once they are completely dry

How to Avoid Spills

Changes in pressure during flights or other travels may cause the medication to leak out of the container. It is recommended to place the applicator in a Ziploc™ or another similar bag to minimize air from getting inside the applicator.

Important Tips

- TICKER® Unlimited dispenses 0.05mL (1/20th of a milliliter) with each CLICK (Assuming there are minimal or no air pockets)
- Gentle Dose® segments are in half-milliliter increments
- TICKER® Unlimited is a one-time-use dispenser and is not refillable. Reusing or refilling it may lead to inaccurate doses
- Gentle Dose® applicator may be reused with the proper care
- Store your medication(s) away from direct sunlight or excessive moisture

Tune into our YouTube channel at [Youtube.com/@Tickerworks](https://www.youtube.com/@Tickerworks) for video DEMOs on the proper use of TICKER® and Gentle Dose® applicator type "Tickerworks Gentle Dose" in the search field.



Please allow your pharmacy at least 48 hours
to prepare your customized medications.