



Seat Color
Black Vinyl

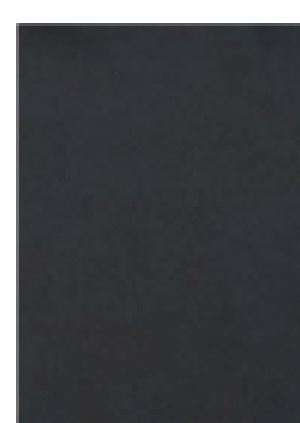

Mossy Oak®
Break Up™ Infinity*

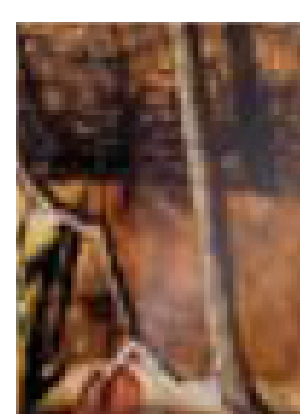

Pink*




Yellow*



White*



Body Colors
Blue
Red

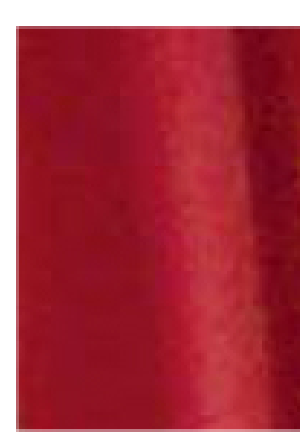

Now available in
more colors!*

*upcharge applies

Ride

iGo Compact

The iGo Compact line expands with the all-new four wheel iGo Compact. With a 300 lb. weight capacity, the four-wheel iGo Compact offers end users more foot room with great indoor maneuverability all with the enhanced stability of a four-wheel scooter. Featuring a Delta tiller with wraparound handles, completely wireless disassembly and two sets of easy-to-change shroud panels in Red and Blue, the four-wheel is the perfect compact scooter for someone looking for a more luxurious ride.



4 wheel specifications

Model	GB146
Medicare Code (pending)	K0800
Weight Capacity	300 lbs.
Drive Wheels	Rear
Maximum Speed	4.5 mph
Operating Range ¹	7.5 miles / 11.0 miles
Ground Clearance	2.5"
Turning Radius	48.5"
Type Batteries	2-12 AH / 2-18 AH
Colors	Red, Blue
Freewheel Mode	Yes
Electronic Speed Control	Yes
Electro-Mechanical Brakes	Yes
Charger	Off-board
Controller	50 Amp Dynamic R-Series
Dimensions	
Length	42"
Width	22"
Height (ground to top of back of seat)	35" -37"
Ground to Top of Seat	20" -22"
Ground to Top of Deck	4.0"
Deck to Top of Seat	15.5" -17.5"
Number of Seat Height Positions	3
Size of Increments	1"
Front Axle to Rear Axle	30"
Weight of Unit: (assembled)	
Front Half	41 lbs.
Rear Half	25 lbs.
Seat w/Arms	21 lbs.
Batteries	21 lbs/pack
Tires:	
Front	8"
Rear	8"
Rear Anti-Tippers	2"
Standard Seat	
	Stadium Style
Back Height - no headrest	16"
Width x Depth	17" x 16"
Color	Black Vinyl
Optional Seats:	
	Stadium Style
Back Height - no headrest	16"
Width x Depth	18" x 16", 20" x 18"
Color	Black Vinyl
Warranty:	
Frame	Lifetime Warranty
Drive Train	2 Year Warranty
Electronics ²	2 Year Warranty

¹ Battery range will vary due to rider weight, drive surface and terrain. Range tested at 200 lbs.

² Electronics warranty excludes batteries.

Wireless disassembly in 4 easy steps!



Step 1: Pop off the seat!



Step 2: Pull out the battery pack!



Step 3: Fold down tiller and lock in place!



Step 4: Pull up on the drivetrain release lever to separate the front and rear frame sections!



Disassembles into 5 pieces!

Features

- Exclusive interlocking teeth for fully connectorless assembly/disassembly
- 300 lb. weight capacity
- Tiller mounted basket