

ScopeValet™ AAMI Level 3 Procedure Gown

Item No. 345EG3

Protection Grade AAMI Level 3: MODERATE FLUID BARRIER PROTECTION

Used for the widest range of surgical procedures, where moderate fluid protection is indicated


ScopeValet™ PROCEDURE GOWNS were designed to meet the specific protection and comfort requirements for endoscopic and other medical procedures. Level 3 Procedure Gown is made of Spunbond-Meltblown-Spunbond (SMMS), a composite fabric designed for the growing needs of the medical industry. This fabric provides increased protection that will not compromise comfort. Level 3 gowns provide waterproof protection for moderate fluid contact.

Spunbond-Meltblown-Spunbond (SMMS) Material Features

- Balanced sheer strength with increased elasticity
- Spunbond comfort with meltblown fluid and micro-organism protection
- Abrasion resistant

Procedure Gown Features


- One size fits all; not individually packaged
- Disposable, one time use
- Latex Free
- Lightweight design and breathable to keep staff cool
- Premium elastic cuffs
- Adjustable Velcro® neck closure to allow for easy on/off and improved usability for workers with mild arthritis or limited mobility
- During removal, the wearer can easily contain contamination by rolling the gown into a ball before discarding it

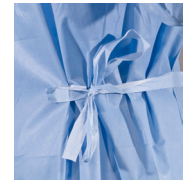

One Size Fits All



Premium Elastic Cuffs

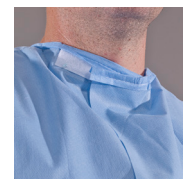

Adjustable Velcro® Neck Closure

Guideline References

AAMI Level 3 Gowns: This classification describes surgical gowns, isolation gowns, other protective apparel, surgical drapes, and drape accessories that demonstrate the ability to resist liquid penetration in two laboratory tests, AATCC 42 (Water resistance: Impact penetration test) and AATCC 127 (Water resistance: Hydrostatic pressure test). For Level 3, the test criterion for AATCC 127 performance has been set at a higher value than for Level 2.