


And then some.



Clock Aneroid Extension

Includes Basket, Coiled Extension Tube, 1 Adult, and 1 Large Adult Inflation System



Fiberoptic Otoscope

Specula Dispenser Extension

Additional Handle Extensions
Each Extension Includes 1 Handle

Add up to 3 Extension Modules

TO CONFIGURE YOUR ADSTATION

Visit our Website at www.adctoday.com or click on the QR code with your smart phone!



Distributor




55 Commerce Drive
Hauppauge, New York 11788
Tel: 1-800-ADC-2670
Fax: 631-273-9659
www.adctoday.com
NEW YORK LONDON TAIPEI NINGBO TOKYO

ADSTATIONLIT
Printed in U.S.A.



Diagnostix™
ADSTATION™
Modular Diagnostic System

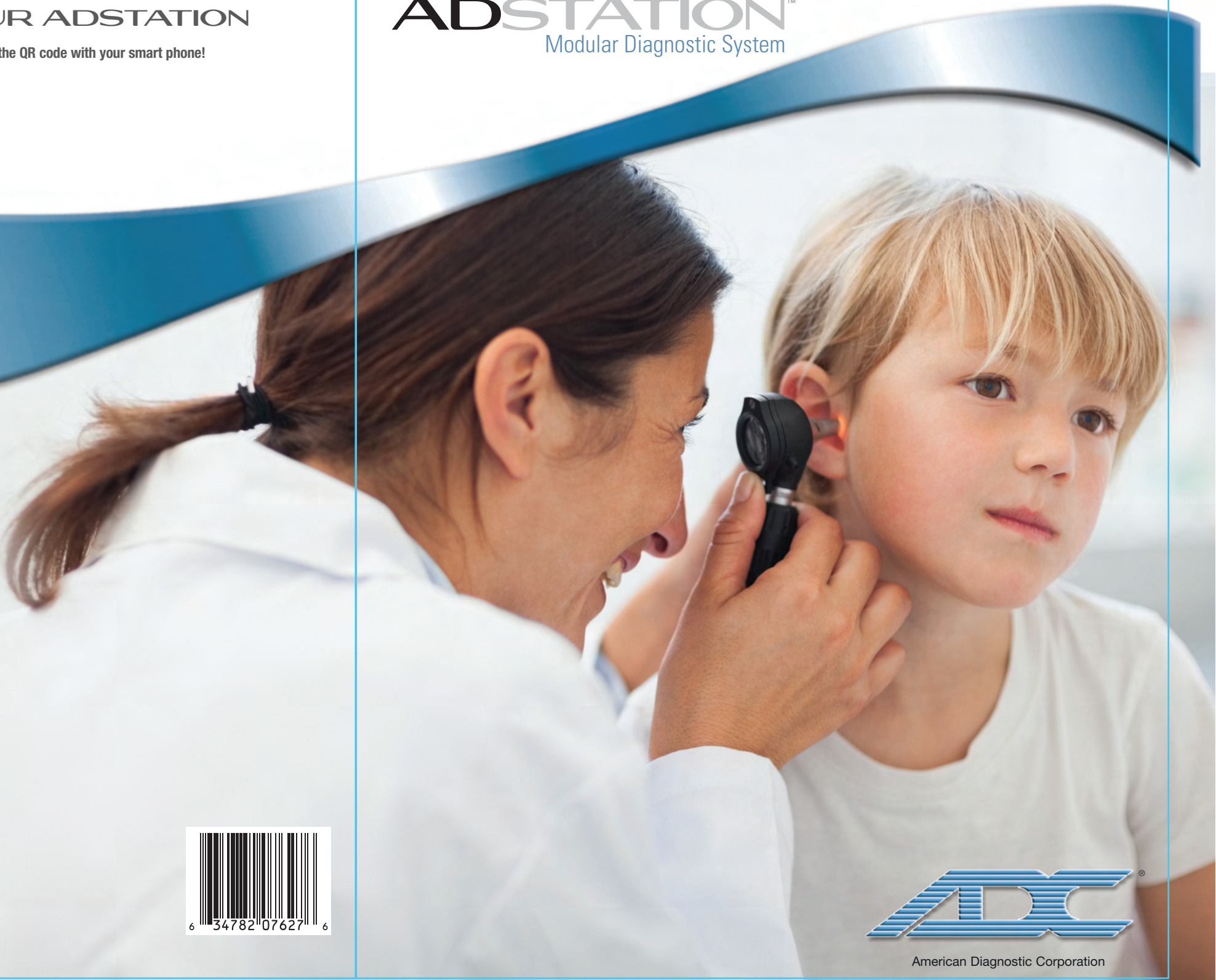

It's everything you've been looking for...

Introducing the ADC® Diagnostix™ Adstation™ Modular Wall system.

The ADC® Diagnostix™ Adstation™ Modular Wall System puts all of the essential diagnostic tools easily within reach. For one exam room or an entire facility, the Adstation™ Modular Wall System is configurable and expandable to include your choice of EENT instruments; manual blood pressure or ADview™ NIBP, temperature and pulse oximetry. Optional EMR connectivity on ADveiw™ available too. The Adstation™ Modular Wall System ensures reliable, cost effective instrumentation, improving workflow, patient care, and your bottom line.

FEATURES

- Expandable and modular
- Start with as little as a single instrument wall transformer and expand with up to 3 more connected extension modules
- Choose from 3.5v Diagnostix Coax or Coax Plus Ophthalmoscope
- 3.5v Diagnostix Fiberoptic Oscope
- 3.5v Throat Illuminator
- 3.5v Dermoscope
- Choose from Halogen/Xenon or AdLED™ lighting technology
- Optional Diagnostix Clock aneroid
- Optional Adview Diagnostic Station with NIBP can be expanded to include Temp and SpO₂ or any combination in up to 8 configurations with optional Bluetooth connectivity
- Optional Wallboard for simplified wall mounting



And then some.



Complete Adview Vital Signs Monitor - BP, Temp, SpO₂

EENT INSTRUMENT HEAD OPTIONS

ADC® Diagnostix™ Instrument heads use a universal bayonette mount for compatibility with other brands of 3.5 volt power supply handles. Anti-theft locking feature is standard. Choose from Halogen/Xenon or advanced AdLED illumination.



Diagnostix™ Coax and Coax Plus Ophthalmoscopes

- High performance Aspherical condenser lens reduces corneal reflection
- Parallel observational and illumination light beams
- Ergonomic design and integral thumb rest improves stability and allows one finger control of all functions
- Magnified/illuminated aperture viewing window
- Coax model with 29 viewing lenses, 6 apertures and 3 filter overlays
- Coax Plus model with 89 viewing lenses, 7 apertures and 3 filter overlays, secondary focusing wheel with zero retention for quick correction of ametropias



Diagnostix™ Oscope:

- Fiberoptic light transmission for an unobstructed view
- Optical glass reduces reflection and provides 3x magnification
- Bi-directional swivel lens allows instrumentation

- Insufflator fitting for pneumatic otoscopy. Compatible with Welch Allyn specula



Diagnostix™ Throat Illuminator:

- PMMA light transmission
- Tongue blade holder facilitates examination of the mouth and throat

- Eject button discards used tongue blade



Diagnostix™ Dermoscope:

- Direct illumination
- Optical grade 10x focusable magnifying lens
- Large textured focusing ring
- Two interchangeable contact plates – one with scale, one without may be disinfected

- Comfortable soft eyepiece may be used with eyeglasses

*Welch Allyn is a Registered Trademark © 2014 Welch Allyn®

AS ALWAYS, ADC PRODUCTS ARE COVERED BY THE MOST COMPREHENSIVE WARRANTIES

- Lifetime warranty on Diagnostix™ clock aneroid, 3 year warranty on inflation system
- 5 year warranty on Diagnostix™ Adstation™ wall transformer and extension modules
- 2 year warranty on Adview™, and instruments
- Lifetime warranty on AdLED lamps