

ELECTROTHERAPY • ULTRASOUND • TRIWAVE LIGHT • THERMOSTIM

DYNATRON
SOLARIS[®]
PLUS

 **Dynatronics**

Cold, Heat, Electrical Stimulation, and Soft-Tissue Mobilization in the palm of your hand.



THERMOSTIM™

The Dynatron ThermoStim Probe™, an optional accessory to our new Solaris Plus, is a soft-tissue mobilization tool that also allows a practitioner to deliver heat or cold therapy in combination with electrical stimulation—all at the same time—significantly reducing the time typically required to administer three separate treatments.

- Temperature Range: 35° - 112° F and reaches target temperature in less than 60 sec.
- Multi-surface head with edges and corners optimizes ThermoStim transfer and tissue mobilization

Power and Flexibility

Solaris Plus delivers 6 Stim waveforms including IFC, Premod, Biphasic, Russian, High Volt, Microcurrent, and Fixed Frequency IFC/Premod. In addition, Solaris Plus offers 3-Frequency Ultrasound and TriWave Light Therapy. **Best of all, Solaris Plus has the power to deliver up to 5 channels of Stim, Ultrasound, and TriWave Light Therapy through both Probe and Pads—all at the same time.**

Exclusive TriWave Light

Deliver any combination of 3-wavelengths of light (red, infrared, and blue) using the Dynatron TriWave™ Light Probe or 2 robust 5" x 7" TriWave Light Pads. Dual Light Pads allow



TriWave Light

for greater coverage of a joint, treatment of 2 separate areas simultaneously, and the ability to treat large areas of the body when placed side-by-side.

#DCP3



#DLP3

#D71CART



DYNATRON
SOLARIS[®]
PLUS



#D719T




The Future is NOW!

The first of its kind, the revolutionary Dynatron ThermoStim™ Probe* combines cold or heat with electrical stimulation and soft-tissue mobilization. Treat a wide variety of conditions, save time, and enhance patient outcomes. ThermoStim, the future of therapeutic modalities.




Target & Target Sweep

By focusing the treatment precisely where it is needed, Dynatronics' unique TARGET feature provides a more effective interferential treatment. Simply glide your finger across the TARGET touch pad to move the center of interference to the site of your patient's pain.

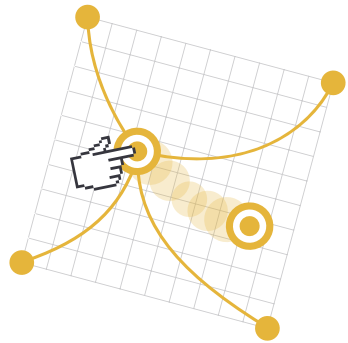

Exclusive 3-Frequency Ultrasound

The only Ultrasound devices in the industry offering 1, 2, and 3 MHz frequencies for the greatest flexibility in depth of treatment.

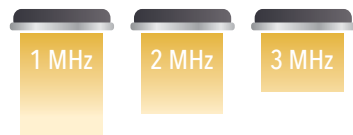

Impressively Simple

Even with nearly unlimited flexibility, setting up a treatment has never been easier. Some treatments can be set up in as few as two key presses. Treatment parameters are highly visible, easy to read, and changed at the touch of a fingertip.

*Patented



#D71CART

Solaris Plus Cart

The Solaris Plus Cart provides ample storage space with 3 large drawers. Constructed with a low center of gravity to prevent tipping, the Solaris Plus Cart is sturdy and aesthetically pleasing.

Dynatron Solaris® Plus Series

The Solaris Plus Series offers the most options in any one device by including Ultrasound, 6 Stim waveforms, and the option of adding TriWave Light Therapy and the ThermoStim Probe. Solaris Plus is covered by a 2-Year Warranty and backed by the best customer service in the industry.

709Plus (#D719T)
5 Channel Combo Unit

708Plus (#D718T)
3 Channel Combo Unit

706Plus (#D716T)
5 Channel Stim Unit

705Plus (#D715T)
3 Channel Stim Unit

Specifications

power requirements _____ 100-240 V~, 50/60 Hz
power consumption _____ 100 W
fuse _____ 250 V, T5AH high break
dimensions _____ 18.5" (w) x 4.0" (h) x 12.0" (d)
47 cm x 10.2 cm x 30.5 cm
weight _____ 8.2 lbs (3.7 Kg)



800 874 6251 | sales@dynatronics.com

www.dynatronics.com

7030 Park Centre Drive Salt Lake City, UT 84121