

HK SURGICAL



aspiration



Pump and
Accessories

Complete Tumescant Solution

HK Aspiration Tubing

The HK Aspiration Tubing is 12 feet in length and has easy grip ends to ensure a secure fit to your machine and handle. This tubing is lightweight, flexible, and extra long for easy maneuvering in the OR. HK Aspiration tubing is compatible with Baxter and Abbott liners. Latex-free.




Order #	Description
AT6-10	10 per case

HK Fine Touch Aspiration Tubing

This tubing is designed specifically to hold luer lock style cannulas without the use of a handle. The built-in luer lock connection allows for easy positioning of the cannula during fine delicate work on areas such as the face and chin. Fine Touch Tubing is 36 inches long and used in conjunction with HK Aspiration tubing. Latex-free.



Order #	Description
ATF-20	20 per case

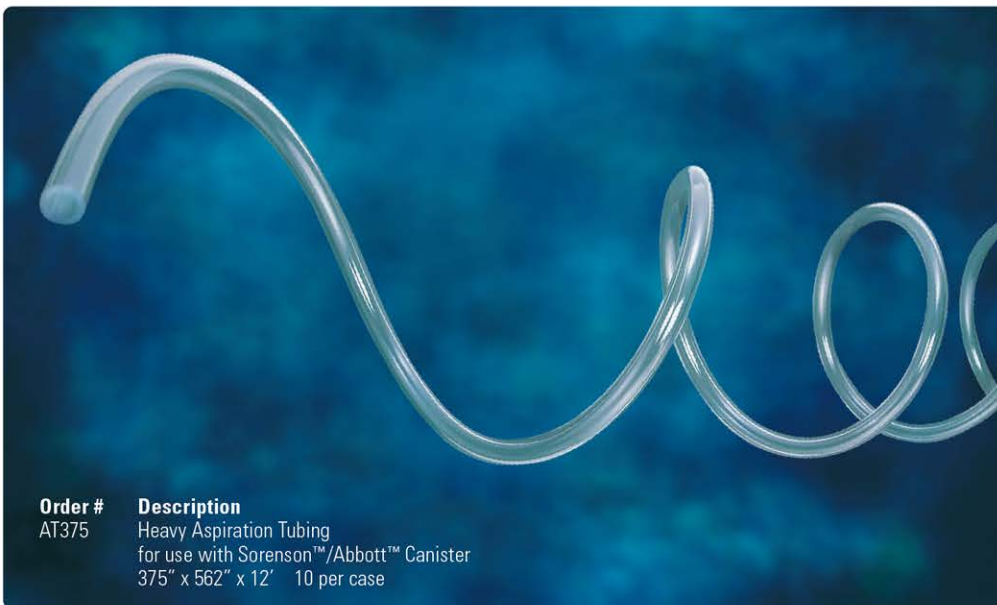

Order #	Description
AT375	Heavy Aspiration Tubing for use with Sorenson™/Abbott™ Canister 375" x 562" x 12' 10 per case

HK Heavy Aspiration Tubing

The HK Heavy Aspiration Tubing is perfect for physicians who prefer a heavier and longer tubing, permitting maximum access around the OR table. HK Heavy Aspiration Tubing is 12 feet in length and compatible with the Sorenson™/Abbott™ canister. Latex-free.



Abbott Liner

The disposable Abbott Liner is used in conjunction with an Abbott Canister. It is soft-shell plastic sealed with a liner lid and can hold up to 1900cc.

Order #	Description
AL-10	10 per case

Abbott Canister

The easy-to-use Abbott Canister is used in conjunction with an Abbott Liner and is made from a durable hard-shell plastic. The Abbott Canister can hold up to 1900cc.

Order #	Description
AC-02	2 per case





HK Finesse Microcannulas

The Finesse Microcannula is designed to minimize the risk of injury to the deep surface of the dermis when liposuction is performed near the skin, making it especially useful during the last stages of liposuction when the dermis is particularly vulnerable to cannula trauma. This microcannula has two in-line oblong slots along one side, and a spatula-like tip.

Order #	Size	
HK10F15	6" (15cm)	10 gauge
HK10F23	9" (23cm)	10 gauge
HK12F10	4" (10cm)	12 gauge
HK12F15	6" (15cm)	12 gauge
HK12F23	9" (23cm)	12 gauge
HK14F10	4" (10cm)	14 gauge
HK14F15	6" (15cm)	14 gauge
HK14F23	9" (23cm)	14 gauge
HK16F5	2" (5cm)	16 gauge
HK16F8	3" (8cm)	16 gauge



The Klein Set

The Klein Set has a combination of Capistrano and Finesse Microcannulas, making it ideal for all body liposuction. Included in the Klein Set are six Capistrano Microcannulas and two Finesse Microcannulas.

Capistrano Microcannulas

Order #	Size	
HK12C15	6" (15cm)	12 gauge
HK12C23	9" (23cm)	12 gauge
HK14C15	6" (15cm)	14 gauge
HK14C23	9" (23cm)	14 gauge
HK16C15	3" (15cm)	16 gauge
HK16C12	4 3/4" (12cm)	16 gauge

Finesse Microcannulas

Order #	Size	
HK12F15	6" (15cm)	12 gauge
HK14F15	6" (15cm)	14 gauge



Order #	Size	
HK10C15	6" (15cm)	10 gauge
HK10C23	9" (23cm)	10 gauge
HK12C15	6" (15cm)	12 gauge
HK12C23	9" (23cm)	12 gauge
HK14C15	6" (15cm)	14 gauge
HK14C23	9" (23cm)	14 gauge
HK16C5	2" (5cm)	16 gauge
HK16C8	3" (8cm)	16 gauge
HK16C12	4 3/4" (12cm)	16 gauge
HK16C15	6" (15cm)	16 gauge
HK18C5	2" (5cm)	18 gauge
HK18C8	3" (8cm)	18 gauge
HK20C6	2 1/2" (6cm)	20 gauge

HK Capistrano Microcannulas

The patented Capistrano Microcannula design was created to maximize liposuction efficiency. Since the apertures are arranged around the entire microcannula circumference, openings are always directed toward the skin.

The Capistrano Microcannula has multiple round apertures arranged in a helical pattern along the distal end of the cannula. It has a single slot located distally to facilitate cleaning and the removal of debris.



The HK Facial Set

The Klein Facial Set has a combination of Capistrano and Finesse Microcannulas, the ideal configuration for chin, cheek, and jowl areas. The set includes two Capistrano Microcannulas, and two Finesse Microcannulas.

Capistrano Microcannulas

Order #	Size	
HK16C8	3" (8cm)	16 gauge
HK20C6	2 1/2" (6cm)	20 gauge

Finesse Microcannulas

Order #	Size	
HK14F10	4" (10cm)	14 gauge
HK16F8	3" (8cm)	16 gauge



Baxter Canister

The Baxter Canister is used in conjunction with a Baxter Liner and is made from a durable hard-shell plastic. The Baxter Canister is easy-to-use and can hold up to 1500cc.

Order #	Description
BC-02	2 per case

Baxter Liner

The disposable Baxter Liner is used in conjunction with a Baxter Canister. It is soft-shell plastic with snap top lid and can hold up to 1500cc.

Order #	Description
BL-20	20 per case



- Vacuum : Adjustable from 0 to 27HG at sea level and standard conditions.
- Flow Rate : 3 cfm free air at sea level and standard conditions. No flow restrictions at intake.
- Dimensions : 19" x 13½" x 16" Approximately 50lbs.
- Limited one year warranty on all parts and labor.
- Includes 1 liner, 1 canister, 2 foot pedals, 1 filter, and 1 piece of aspiration tubing.



HK Aspirator

The HK Surgical Aspirator offers superior high-speed powered suction performance. It can accommodate up to four Baxter Canisters or two Abbott™ / Sorenson™ Canisters, eliminating time-consuming canister changes. The HK Surgical Aspirator is designed for powerful performance with minimal noise.

Order #	Description
AP-III	Aspiration Pump
AP230-III	Aspiration Pump 230 volt

HK Aspiration Filter

The HK Aspiration Filter filters the air and extends the life of your aspirator. It can prevent the accidental intake of fat into the aspirator motor.



Order #	Description
AF-6	6 per case

HK Adit / Skin Punches

A micro-adit used in tumescent surgery is a small circular hole made by a tiny skin biopsy punch. Adits facilitate and promote the open drainage of residual blood-tinged anesthetic solution associated with tumescent surgeries. Skin punches also allow the use of microcannulas with virtually no epidermal friction. The punches are available in 3 sizes 1mm, 1.5 mm and 2mm.

Order #	Description
SP-1.0	1.0mm Skin Punches 50 per case
SP-1.5	1.5mm Skin Punches 50 per case
SP-2.0	2.0mm Skin Punches 50 per case



The Klein Cart

The Klein Cart is specially designed to accommodate and organize your entire system. The top shelf is ideal for holding surgical instruments, while keeping them in a sterile field, during a procedure. The second shelf is for the Klein Infiltration Pump, while the lower area accommodates the HK Aspirator and features "S" hooks for hanging IV bags. Removable trays and locking wheels add stability and mobility.

Order #	Description
KC	Klein Cart



Aspiration Handles

All of the HK Surgical aspiration handles have a luer lock fitting (so they may be used with any luer lock cannula) and can be ordered with or without a hole. The handles are available in two different materials, stainless steel or aluminum; holes are available on either metal. The base of the handle can accommodate all styles and sizes of aspiration tubing.



Order #	Description
AHL	Aluminum Handle Large
AHLH	Aluminum Handle Large with hole
AHS	Aluminum Handle Small
AHSH	Aluminum Handle Small with hole
HKHAND	Klein Aspiration Handle
HKHANDH	Klein Aspiration Handle with hole