

UltraCare® XT

Leading bed evolution.



Joerns® UltraCare® XT

The standard UltraCare XT has the same great features you enjoyed with the UltraCare bed plus newly-designed controls, faster positioning speed, AdvancedCare and IPX4 fluid ingress protection rating to name a few advantages. Then, upgrade with various options and accessories, including UltraWide and UltraLock®. One bed. A total solution.

Features

- 1 Toroidal architecture** engineering of the bed's steel frame is up to 30 times stronger than angle iron
- 2 Therapeutic grid surface** offers improved air flow and ultimate support for the sleep surface
- 3 Integrated pendant connections** on both sides of mattress support platform
- 4 Four-section mattress support platform** offers full range of comfort and therapeutic positions, which may reduce skin shear
- 5 Attractive bed panel options** available in a wide array of designs and finishes
- 6 Staff control** allows caregivers flexibility to select which features may be adjustable by the resident
- 7 EasyLift Leg Lift** provides quick positioning for improved circulation
- 8 Revolutionary arched legs** deliver a super-strong, wide stance for equipment clearance
- 9 Four DC motors** supply amazingly quiet and smooth operation
- 10 7" to 30" travel range** offers safety, whether resting at an ultra low height or providing care at the higher end
- 11 Assist devices** have superior functionality and perfect grip height
- 12 UltraLock** is available at the foot end or both the head and foot to maximize security or mobility at any height (option)
- 13 AdvancedCare** features gravity-assist for easier and safer repositioning. Utilize Comfort Chair functionality for comfort, relaxation, improved circulation, and more ergonomic care by bringing the resident closer to the caregiver (not shown).




Simplified Controls

Ergonomically designed pendant meets geriatric needs featuring soft touch edges, finger grips and optional enhanced back lighting, raised buttons and illuminated lock-out status.

Options and Accessories

The UltraCare XT can be easily transformed to fit any configuration or room style. Upgrade to suit your facility's needs. Choose options such as SafeLight, UltraLock and UltraWide or experience other great features to customize your UltraCare XT bed.

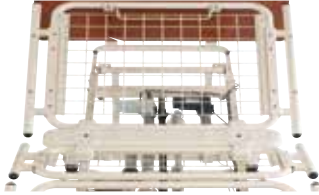

ULTRAWIDE

Transform your UltraCare XT to the widest 500lb (227kg) safe working load bed on the market. The extra wide width makes for safer repositioning and a more comfortable night's sleep. Choose from 39" or 42" extension kits.



ULTRALOCK

Designed for mobility or security at any height, while allowing for a 7" low height. UltraLock provides quick, easy bed movement.



MOBILITY STATUS INDICATOR

Visually alerts your caregivers with a flashing green light that the bed mobility is in an unlocked position. Easily attaches to each caster, adding a built-in safety reminder to each bed.



PENDANT HOLDER

Smooth design easily affixes to Joerns assist devices making pendant stowage and retrieval simple. Residents can even operate the pendant while sturdily stowed out of the way.



SAFELIGHT

Available with three light settings and ability to install without tools, SafeLight illuminates the underside of the bed on either side while providing light during the evening for individual ingress/egress or safer transfers.



4" AND 8" LENGTH EXTENDER

Extend the length of your UltraCare XT mattress surface by 4" for 35" width frames or 8" for 42" width frames to accommodate taller, larger residents.

Assists



ASSIST DEVICE

Features a low profile design, avoiding feeling of restriction or restraint for the resident. Mounted at the perfect ergonomic height to assist residents in and out of bed.

Soft-touch handle provides a sturdy handhold for improved resident grip.



DELUXE ASSIST HANDLE

New assist handle features an ergonomic and a sturdy soft touch handhold.

The assist is mounted at an ergonomic height to best assist residents in and out of bed, and has an integrated pendant holder on each side.

Technical Specifications

| | |
|--|------------------|
| Safe Working Load | 500 lbs / 227 kg |
| Under Bed Clearance | 23.7" / 603 mm |
| Product Weight | 220 lbs / 99.7kg |
| Max Back Angle to Horizontal | 65° |
| Max Knee Angle to Horizontal | 35° |
| Min Height From Floor (Surface Platform) | 7" / 178 mm |
| Max Height From Floor (Surface Platform) | 30" / 762 mm |
| Width (minimum) | 35" / 889 mm |
| Width (maximum) | 42" / 1067 mm |
| Length (minimum) | 76" / 1930 mm |
| Length (maximum) | 88" / 2235 mm |

Electrical Specifications

| | |
|----------------------------|-----------------------------|
| Mains Voltage | 110-220 VAC |
| Frequency | 60 Hz |
| Input Current | 4.0 Amps |
| Control Box Output Voltage | 24V DC |
| Electrical Protection | Class I or II, Type B |
| Duty Cycle | 10% (2 min on / 18 min off) |
| Circuit Protection | Auto reset current sensor |



Frame Construction

- All major structural components, including support surface platform, are designed with steel tubing in combination with patented toroidal construction to provide the utmost in structural stability

Additional Information

- All exposed steel surfaces are coated with epoxy powder paint
- Fasteners, screws, rivets and other non-painted parts are corrosion resistant
- All control switches are of momentary contact design
- Head panel is easily lifted off the bed to allow for CPR in emergency situations

U.S. Patents

- 6473922, 6474869, 6519794, D436,474, 437,143, D445,614
Other Patents Pending

Finish

Joerns UltraCare XT comes in a soft tone finish.

We reserve the right to make dimension or specification changes, or deletions consistent with design trends, product improvements, or for reasons beyond our control.

Color options have been reproduced as faithfully as printing processes allow.

Certifications — ETL Listed

- USA: Conforms to UL 60601-1, 1st ED., 2003
- Canada: Certified to CAN/CSA C22.2 No. 60601.1-M90, CAN/CSA STD C22.2 NO. 60601-2-52-03
- International: Conforms to IEC 60601-2-52, Third ED., 2009, IEC 60601-2-52, Amendment 1.1, 2016
- European: EN60601-1, EN60601-2-52

Environmental Conditions:

Ambient Temperature: 50°F to 104°F / 10°C to 40°C
Relative Humidity: 30% to 75% (non-condensing)

Classification:

Conforms to: IEC 60601-2-52

