



A leader in infection prevention
since 1925

**HIGH THERMAL DISINFECTION WASHER
TIVA8-MD | 165 LITER CHAMBER**

Medical



Quality Products ♦ Leading Edge Technology ♦ Reliable Services ♦ Professional Excellence

Welcome to the TIVA8-MD

High Disinfection Thermal Washer

The ultimate automated solution for washing, rinsing, disinfecting, and drying your complex instruments prior to sterilization. Designed to simplify your work by eliminating manual washing and reducing your manpower. Ensuring the safety of your patients by providing consistent quality washing cycles at all time and ensuring the safety of your staff by eliminating risks of exposure to potentially harmful aerosols.



Ambulatory surgery centers



Hospitals CSSD




Physician's Office



Vet hospitals/clinics

TIVA8-MD Washer Disinfector Package with Drying

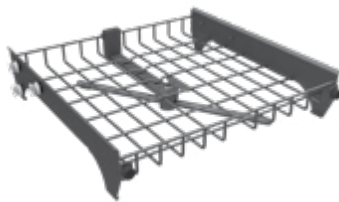

Upper Rack with spray arm



Lower Rack basic

Features included:

- Double-layered tempered glass door
- Three water connections: Cold + Warm + Deionized water softener with protection filter
- Spray Arm Monitoring
- Multi Color Chamber LED Light
- Upper Rack with spray arm
- Lower Rack basic

Comes in 2 different power configurations:

- 3-Phase 208-240V without Neutral, 60Hz
- Single Phase 208-240V without Neutral, 60Hz

Optional Medical Package:

TIVA8 Medical Package P/N:
ATS238-0524



4 X Din Baskets (19" x 10" x 2")

Improved Productivity and Efficiency



Premium Performance

- Instruments are completely dried thanks to powerful electric blowers and heaters with a HEPA H14 filter
- Three water connections for flexible operation and for optimal washing needs (cold & hot with softeners, demineralized)
- Conductivity measurements to check rinse water quality before disinfection
- Up to 40 washing programs
- AISI 316L Stainless Steel



Intuitive User Experience

- A full control and all-in-one operation from a single touch panel with color display
- Spray Arm monitoring for continuous washing validation
- Graphical report of washing phases including time and temperature
- Multi color chamber LED light (White, Green, Red) operates during loading, unloading and running the program
- Double-layered tempered glass door for visual monitoring



Enhanced Traceability

- Comply with regulation by saving washing cycles data
- USB and RS232 connectivity
- Printer for data documentation (optional)




Eco-Friendly

- 2 dosing pumps with flowmeter for precise consumption of acid and alkaline detergents
- Integrated water softener for reducing decontamination
- Built-in steam condenser for excess water and steam

TIVA8-MD technical data

23.6 inch narrow design

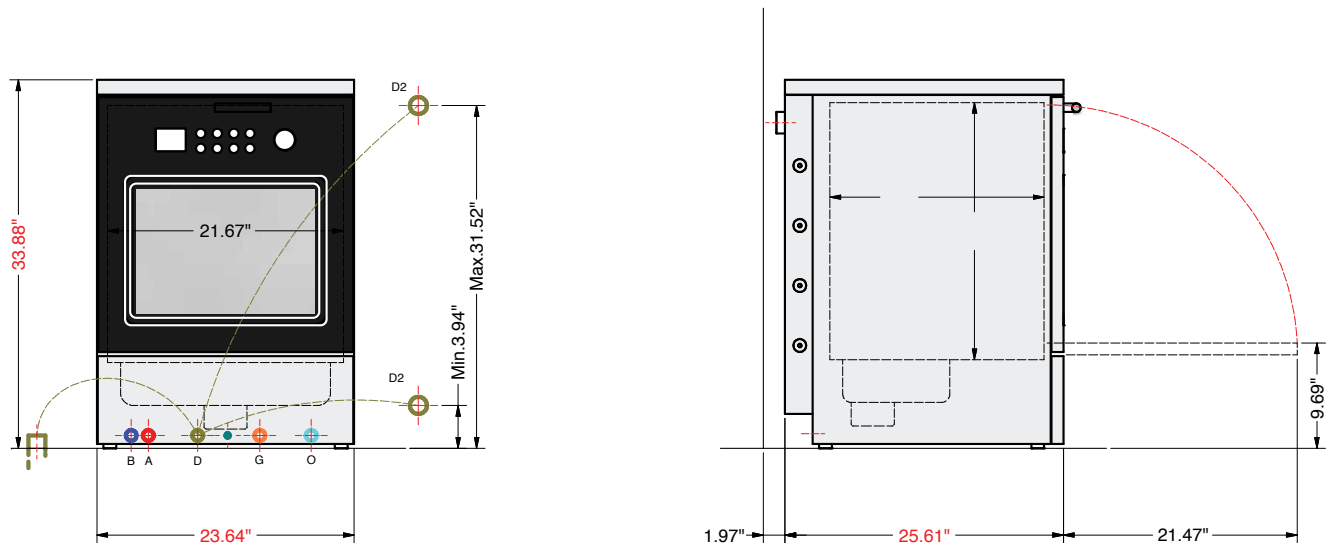

A Hot water connection-3/4" (tube 1/2" g)

B Cold water connection-3/4" (tube 1/2" g)

D2 Drain dn 40 - hose connection 0.984" (at wall)

G Electrical connection

O Demineralized water connection-1/2" g

E2 Drain dn 40 - hose connection 25mm (at wall)

ALL YOUR INFECTION CONTROL NEEDS

Let us introduce you to our complementary family of sterile processing products to ensure a comprehensive process from start to end



T-Clean TIVA Detergent:

Detergent for reprocessing thermostable and thermolabile instruments



T-Clean TIVA Neutralizer:

Acidic neutralizing agent and detergent for the automated reprocessing of instruments and laboratory glassware



T-Clean Prep-Foam:

Multi-Tiered Enzymatic Foam Spray. Helps To Prevent The Adhesion Of Bioburden To The Surfaces Of All Surgical Instruments And Scopes



T-Edge



Water Distiller



Ultrasonic Cleaner



Cleaning Indicator

LTR330-0109/WAS/TIVA8MED/USA/V-1/02.22

Medical

Distributed by:

Tuttnauer
Innovation · Legacy · Partnership

Tuttnauer USA Co.

25 Power Drive,
Hauppauge, NY 11788
Tel: +800 624 5836, +631 737 4850
Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com



Learn from our experts
Join our blog tuttnauerusa.com/blog