



- **NovaSeal®** for safer separation to ensure sterility

- Polypropylene Luer Locks connect small diameter tubing with secure, leak-proof connections.



- Polypropylene Plug is used to keep inner tubing free of contaminants before use.

- CPC MPC connector design allows for quick sterile connections and the ability to transfer media easily.



- CPC AseptiQuik® Connectors provide quick and easy sterile connections for small-flow applications. The “FLIP-CLICK-PULL” design enables users to quickly transfer small volumes of media easily with less risk of operator error.


MADE IN

Foxx Life Sciences EZBio® Made to Order “MTO” Single Use Sterile Media Bottle Assemblies with Tubing, and Vent Filters. Bottles available in Polyethylene Terephthalate Glycol-modified (PETG), or Polycarbonate (PC). **EZBio®** “MTO” Single Media Bottle Use Assemblies, allow the customer to select individual components suited for their project needs. Selectable components include Volume and Material of Bottle, Type and Brand of Fittings with Male and Female connectors and Brand of TPE Tubing’s with NovaSeal.

All **EZBio®** “MTO” Single Use Media Bottle Assembly, include, Bottle, Cap, Cap Adapter (Insert), Extra Close cap for storage, TPE Tubing ID 1/8in (3.2 mm) X OD 1/4 in (6.4 mm), Fittings with Male and Female connectors, a Dip tube, for an accurate draw, a 0.22 μ m Vent filter, Ratchet Style Pinch Clamp, and an option to include **NovaSeal®** for quickest aseptic tube-separation. Each finished is double bagged and Gamma Sterilized in the US.

Separating one part of a disposable assembly from the other is often required during equipment operation in the pharmaceutical industry. When it comes to biopharmaceuticals, the separation and isolation of compartments can cause cross-contamination. The **NovaSeal®** is a metallic pinch pipe on the exterior of the tubing. It provides a sterile means for separating a disposable assembly by crimping or cutting a plastic tool at a given location. It ensures that sterility is maintained and the fluid is secured and protected from outside.



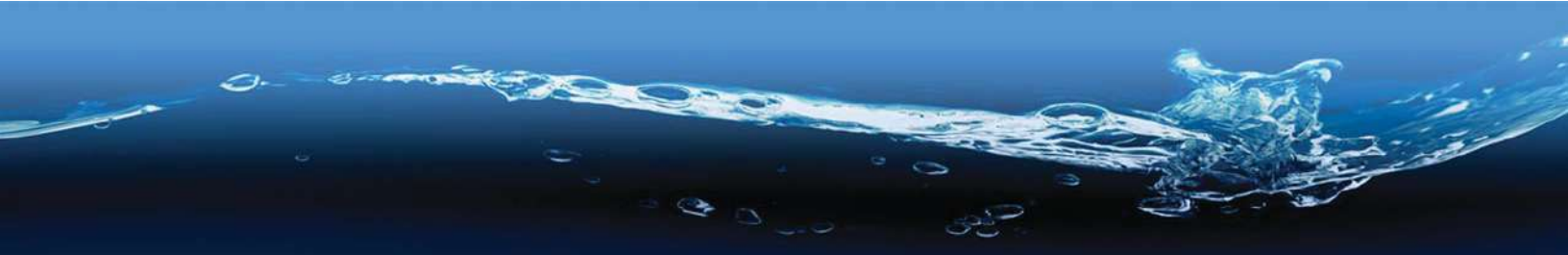

BUILD YOUR MTO

43X-XX0X-FLS

Product Line

43: EZBio® Single-Use Kits

Bottle Material

0: Cap Assembly (No Bottles)
 7: PETG
 8: PC

Tubing Brands

8: EZBio® Tubing w/ NovaSeal
 9: C-Flex® Tubing w/ NovaSeal
 A: AdvantaFlex® Tubing w/ NovaSeal
 B: Imaweld® Tubing w/ NovaSeal
 C: Flexelene™ Tubing w/ NovaSeal
 D: PureWeld® Tubing w/ NovaSeal

Bottle Volume

2: 60mL
 3: 125mL
 4: 250mL
 5: 500mL
 6: 1000mL
 7: 2000mL

Fittings

0: Plug
 1: MPC CPC
 2: Luer
 3: AseptiQuick® S

