## **EZBio® MTO Cap Assemblies**

www.foxxlifesciences.com



- · Four different connectors available
- Polypropylene Luer Locks connect small diameter tubing with secure, leak-proof connections.





 Polypropylene Plug is used to keep inner tubing free of contaminants before use.

 CPC MPC connector design allows for quick sterile connections and the ability to transer media easily.





· CPC AseptiQuik® Connectors provide quick and easy sterile connections for small-flow applications. The "FLIP-CLICK-PULL" design enables users to quickly transfer small volumes of media easily with less risk



of operator error.



EZBio® Cap Assemblies are made entirely with USP Class VI components and Gamma Irradiated in the United States. Assemblies take advantage of Foxx's patented VersaCap® technol-ogy which means there is no machining necessary to produce the product. Perfect for single-use applications and to phase out costly stainless steel apparatuses expensive to clean and validate.

All EZBio® MTO Cap Assemblies include a polypropylene cap, cap adapter, TPE tubing ID 1/8in (3.2 mm) X OD 1/4 in (6.4 mm), fittings with Male and Female connectors, a dip tube, for an accurate draw, a 0.22  $\mu$ m vent filter and a ratchet style pinch clamp. Each finished is double bagged and Gamma Sterilized in the United States.











## **EZBio® MTO Cap Assemblies**

www.foxxlifesciences.com



## YOUR EZBIO® MTO CAP ASSEMBLY. YOUR CHOICE.




## **Ordering Instructions**

Foxx Life Sciences has made it easy for its customers to have full customization when ordering EZBio® MTO Cap Assemblies on **ShopFLS.com**.

PICK YOUR TUBING

PICK YOUR FITTING

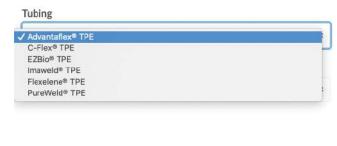

Tubing Fitting

✓ CPC AseptiQuik®

Luer

CPC MPC
Plug















