

- **VersaCap®** technology provides the flexibility of interchangeable top-connectors
- **VersaCap®** technology is offered in 50mm, 80mm, and 120mm industry standard sizes
- Quick-connect adapters available for easy port changes
- Leak-proof cap with replaceable gasket
- **VersaCap®** design enables stacking
- Available in a variety of configurations
- **VersaCap®**, **VersaBarb®**, and Spigots are autoclavable
- **VersaBarb®** offers a convenient solution for bottom-connections and can fit a variety of hose sizes
- Spigot offers two flow rates

Foxx Life Sciences has enhanced the features and benefits of the **EZgrip™** carboys and Heavy Duty Vacuum Bottles by offering the **VersaCap®**, **VersaBarb®**, and state-of-the-art spigot technology. These three components expand the wide array of applications supported by the **EZgrip™** carboy.

Patented **VersaCap®** is leak-proof, has replaceable platinum-cured silicone gaskets, and provides the ultimate in top-connection flexibility. Open **VersaCap®** is compatible with a variety of adapter inserts, enabling users to quickly change the top connections of each carboy or bottle to match their application needs. **VersaCap®** offers large handles with finger grips that allow for easy tightening and opening. They are also stackable and are offered in 50mm, 80mm (compatible with most carboys), and 120mm wide-mouth versions that fit exclusively on Foxx Life Sciences' 20L, 40L, and 75L **EZfil™** carboys.

Quick-connect hose barb options include quick-connect port combinations of two, three, and four, with  $\frac{1}{4}$ " and  $\frac{3}{8}$ " hose barb sizes, as well as molded-in hose barb options, available with  $\frac{1}{8}$ ",  $\frac{1}{4}$ ",  $\frac{5}{16}$ ", and  $\frac{1}{2}$ " sizes. The Clearview® adapter is also available for users who need a crystal-clear window into the carboy for visualizing solution homogenization.


For bottom tubulation options, **VersaBarb®** provides flexibility in drainage, pump connection, filling bell, and automated equipment. **VersaBarb®** is available in  $\frac{1}{4}$ ",  $\frac{3}{8}$ ",  $\frac{1}{2}$ ", and  $\frac{3}{4}$ " hose barb sizes. A  $\frac{3}{4}$ " sanitary fitting **VersaBarb®** is also available to enable the **EZgrip™** carboy to function as a supply reservoir to a larger system.


Foxx Life Sciences' **EZgrip™** carboy spigots have large levers that control two flow rates, one for continuous, locked flow and the other for a faster, manual flow. The spigot spout features  $\frac{5}{8}$ " and  $\frac{3}{4}$ " hose barbs, and a standard internal luer tape to fit a male barb adapter. The spigot body has internal threads with a hard stop to prevent over-tightening which may damage the seal.


All **VersaCap®**, adapters, **VersaBarb®**, and spigots are made from polypropylene, which may be autoclaved at 15 PSI and 121°C (250°F) for 20 minutes.






**VersaCap® 80mm**


<b>205-1003-RLS</b>				
	Cap, Closed			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	n/a	n/a	2


<b>205-1004-RLS</b>				
	Cap, Open with closed adapter			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	n/a	n/a	2


<b>205-1007-RLS</b>				
	Cap, Open top with adapter, 1/4" HB x 2 and Vent			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	n/a	1


<b>205-1008-RLS</b>				
	Cap, Open top with adapter, 1/2" HB x 2 and Vent			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/2"	n/a	1


<b>205-3001-RLS</b>				
	Adapter, Closed			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	n/a	n/a	1


<b>205-3002-RLS</b>				
	Adapter, 1/4" HB x 2			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	2	1


<b>205-3003-RLS</b>				
	Adapter, 3/8" HB x 2			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	3/8"	2	1


<b>205-3004-RLS</b>				
	Adapter, 1/4" HB x 3			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	3	1


<b>205-3005-RLS</b>				
	Adapter, 3/8" HB x 3			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	3/8"	3	1

<b>205-3006-RLS</b>				
	Adapter, 1/4" HB x 4			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	4	1

<b>205-3007-RLS</b>				
	Adapter, 3/8" HB x 4			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	3/8"	4	1

<b>205-5001-RLS</b>				
	Adapter, 1/4" HB x 2 and Vent			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/4"	n/a	1

<b>205-5002-RLS</b>				
	Adapter, 1/2" HB x 2 and Vent			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	1/2"	n/a	1

<b>208-3001-RLS</b>				
	Adapter, Closed Clearview®			
	Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
	80	n/a	n/a	1

HB = Hose Barb


**VersaCap® 50mm**

**205-0313-RLS**

Cap, Closed

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	n/a	n/a	2


**205-0345-RLS**

Cap, Vented top with 0.2 µm Hydrophobic Membrane

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	n/a	n/a	1


**205-0314-RLS**

Cap, Open with closed adapter

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	n/a	n/a	2


**205-3331-RLS**

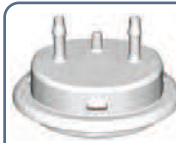
Adapter, Closed

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	n/a	n/a	1


**205-0317-RLS**

Cap, Open top with adapter, 1/8" HB x 2 and vent

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	1/8"	n/a	1


**205-5331-RLS**

Adapter, 1/8" HB x 2 and vent

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	1/8"	n/a	1


**205-0318-RLS**

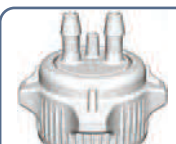
Cap, Open top with adapter, 1/4" HB x 2 and vent

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	1/4"	n/a	1


**205-5332-RLS**

Adapter, 1/4" HB x 2 and vent

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	1/4"	n/a	1


**205-0319-RLS**

Cap, Open top with adapter, 5/16" HB x 2 and vent

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	5/16"	n/a	1


**205-5333-RLS**

Adapter, 5/16" HB x 2 and vent

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
50	5/16"	n/a	1

HB = Hose Barb




**VersaCap® 120mm**

**205-2003-RLS**

Cap, Closed

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	n/a	n/a	2


**205-2004-RLS**

Cap, Open with closed adapter

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	n/a	n/a	2


**205-4001-RLS**

Adapter, Closed

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	n/a	n/a	1


**205-4002-RLS**

Adapter, 1/4" HB x 2

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	1/4"	2	1


**205-4003-RLS**

Adapter, 3/8" HB x 2

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	3/8"	2	1


**205-4004-RLS**

Adapter, 1/4" HB x 3

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	1/4"	3	1


**205-4005-RLS**

Adapter, 3/8" HB x 3

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	3/8"	3	1


**205-4006-RLS**

Adapter, 1/4" HB x 4

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	1/4"	4	1


**205-4007-RLS**

Adapter, 3/8" HB x 4

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	3/8"	4	1


**205-4008-RLS**

Adapter, EZ Top® adapter

Cap Size	Hose Barb Size	No. of Quick Connects	Qty/pk
120	n/a	n/a	1




**VersaBarb® and Spigot Technology**

**225-4001-RLS**

VersaBarb®, Spigot fitting, closed cap

Hose Barb Size	Thread Diameter	Qty/pk
n/a	1 1/8"	1


**225-4005-RLS**

VersaBarb®, Spigot fitting, with 3/4" HB

Hose Barb Size	Thread Diameter	Qty/pk
3/4"	1 1/8"	1


**225-4002-RLS**

VersaBarb®, Spigot fitting, with 1/4" HB

Hose Barb Size	Thread Diameter	Qty/pk
1/4"	1 1/8"	1


**225-4006-RLS**

VersaBarb®, Spigot fitting, with 1/2" Sanitary Connector

Hose Barb Size	Thread Diameter	Qty/pk
n/a	1 1/8"	1


**225-4003-RLS**

VersaBarb®, Spigot fitting, with 3/8" HB

Hose Barb Size	Thread Diameter	Qty/pk
3/8"	1 1/8"	1


**225-1001-RLS**

Spigot, 2 Pack

Hose Barb Size	Thread Diameter	Qty/pk
n/a	1 1/8"	2


**225-4004-RLS**

VersaBarb®, Spigot fitting, with 1/2" HB

Hose Barb Size	Thread Diameter	Qty/pk
1/2"	1 1/8"	1

HB = Hose Barb



We recommend testing for your specific application before use.