

Hardwire Pairing Kit

Quickstart Guide

Hardwire Pairing Kit Multiple Device Setup

Half Max, Duo, Max, Quad, Elite

To pair your Multiple Device Setup using Pairing Cable Kit follow the steps listed below.

Reset All Devices

Step 1

After applying power to the device using the main AC power cord, reset device by pressing and holding the Power and Light Mode buttons together for 4 seconds.

Step 2


Press the Power button twice or until the display reads FOLL.

Step 3

Press and hold the Power button for 4 seconds to confirm.

Unplug the power cord and leave the device in this state. The next time the device receives power it will search for the Lead device.

Repeat steps 1-3 for each device before continuing.

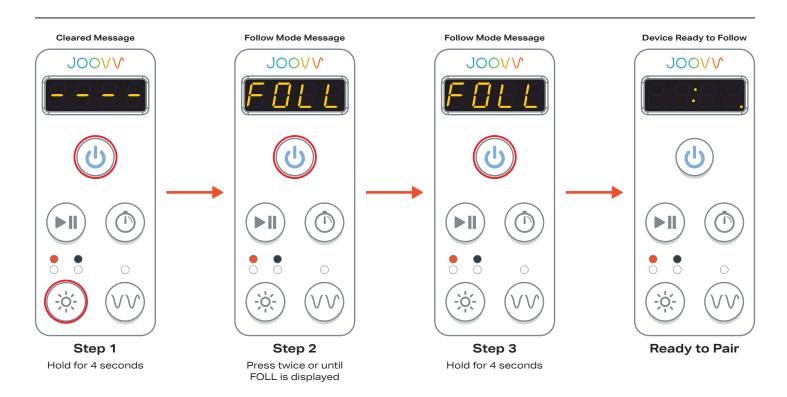

Pairing Cable Kit steps for Multiple Device Setup

Assemble Joovy

Step 4

Assemble and mount devices then connect AC patch cables, main AC power cord, and hardwire pairing cables following the Wire Diagram section on pages 29-32 of the User Manual.

Set Lead Device

Step 5

Locate the Lead device control panel (Example Joovv Quad shown). **Press and hold the Power and Light Mode buttons for 4 seconds** to clear the device mode.

Step 6


Press the Power button once or until the display reads LEAD.


Step 7

Press and hold the Power button for 4 seconds to confirm. Devices will automatically pair and the Lead control panel will display the 10:00 minute session timer.

The Lead device is now the controllable device for the configuration.

Great job! You're now ready to start Joovvin'!

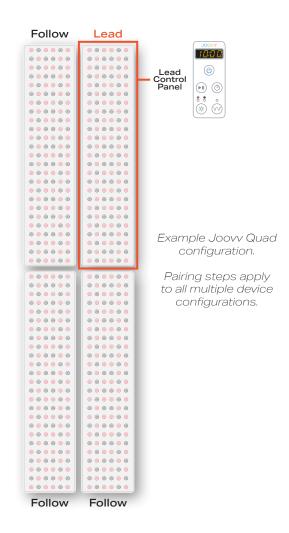


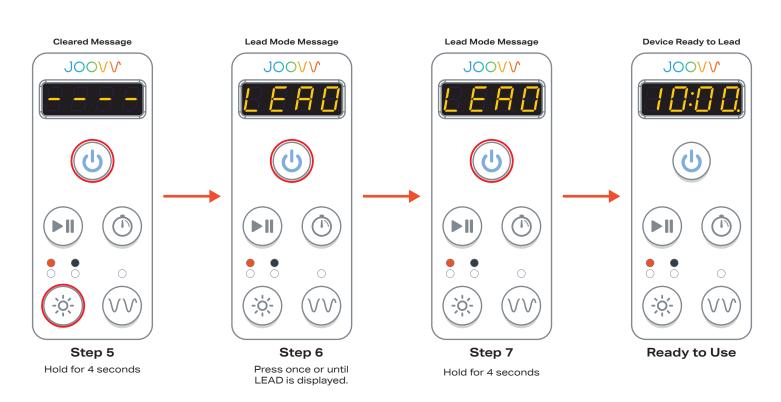


Device Configurations

with pairing cables displayed

Joovy Half Max 3.0

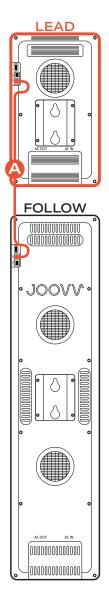

Joovv Duo 3.0

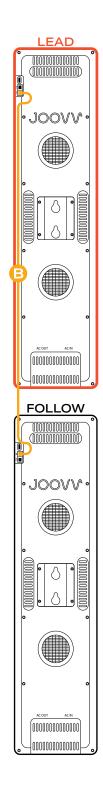


Device Configurations

with pairing cables displayed


Joovv Duo 3.0(horizontal)

Joovy Max 3.0

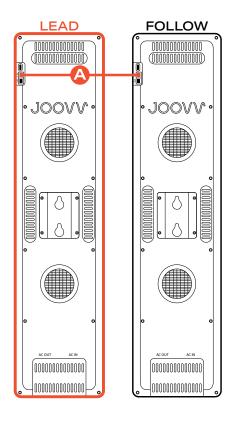




The Horizontal Duo configuration is an alternative setup option and requires an additional accessory kit. Please reach out to help@joovv.com for more information.

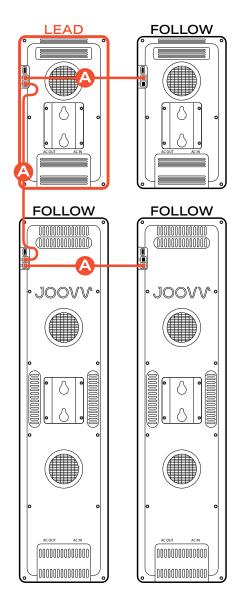


Device Configurations

with pairing cables displayed

Joovy Quad 3.0

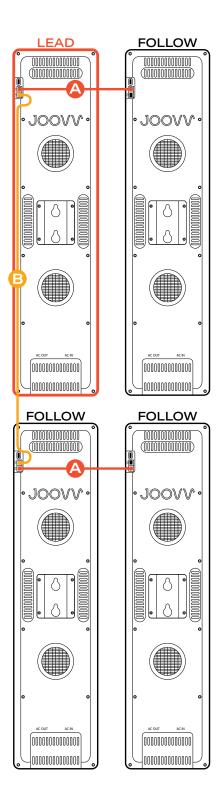

Joovy Elite 3.0

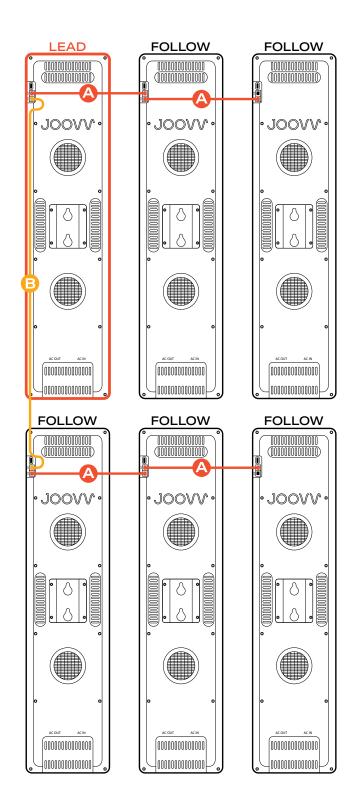

Revision History		
Date	Revision #	Description of Changes
December 2020	JV-MKT-DOC-3523 (Ver.0)	Ver.0 Release to market
March 2021	JV-MKT-DOC-3523 (Ver.1)	Ver.1 Updated Page 2 Illustration
April 2021	JV-MKT-DOC-3523 (Ver.2)	Updated text/images

