

Help Alert




Quick Start Guide


Table of Contents

1. Front Layout
2. Back Layout
3. Side Layout
4. What Do the Lights Mean?
5. Connecting Magnetic Cable
6. Connecting Charging Base
7. Plug into AC outlet
8. Charging the Device
9. Turning Device "ON"
10. Placing a Call
11. Using the Speakerphone
12. Disconnecting a Call
13. Turning Device "OFF"
14. Turn Voice Prompts On/Off

Front Layout

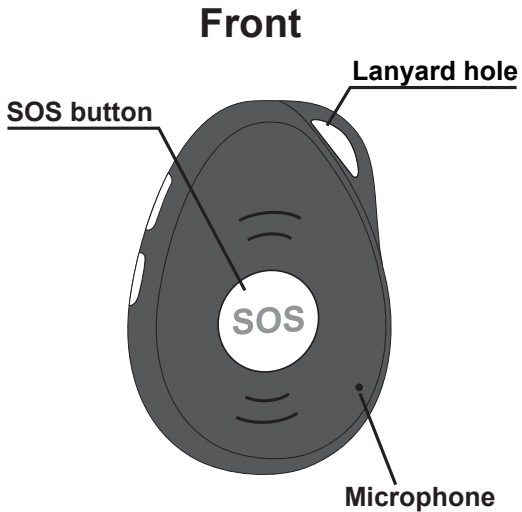


Back Layout

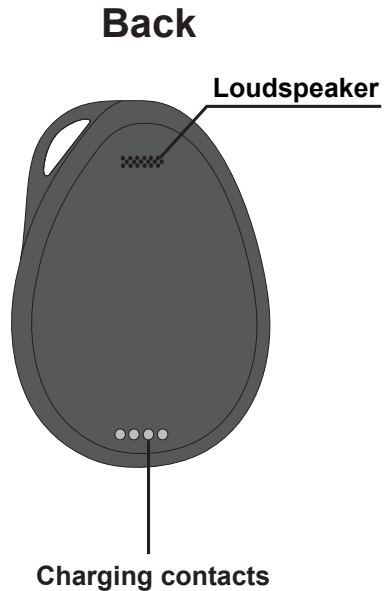

Side Layout

Left Side

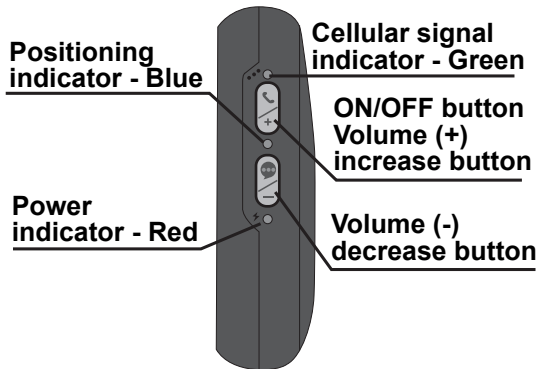

What Do the Lights Mean

Cellular signal indicator - - Green

Green	Light shows a single flash rapidly every 3 seconds	Light shows a double flash rapidly every 3 seconds
Means	The device has a stable Cellular signal	The device is registered to the Cellular network

Positioning indicator - - Blue

Blue	Light shows a single flash rapidly every 3 seconds	
Means	The device has no latest location fix	


Blue	Light shows a double flash rapidly every 3 seconds	Light Off
Means	The device has latest location fix	The device is not fixing the latest location

Power indicator - - Red


Red	Red ON (solid)	Red shows a double flash rapidly every 3 seconds
Means	Device has been fully charged	BLE connected

Red	Red Blinking Quickly	Red Off or blinking slowly
Means	Battery power is lower than 20%	The device is charging

Connecting Magnetic Cable

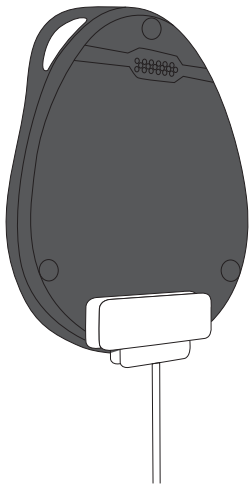


Connecting Charging Base

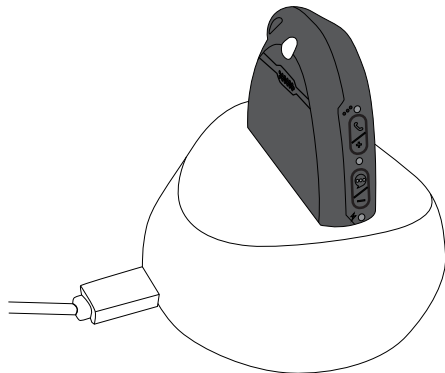


If you are using the optional charging stand, plug the USB charging cable into the charging base.

Plug into AC Outlet

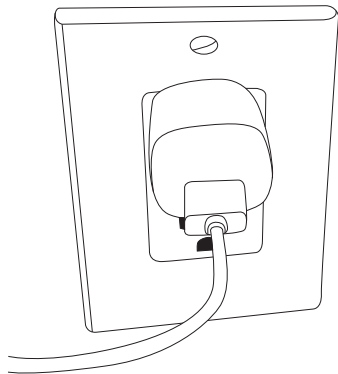


Plug the charging cable or charging stand into the wall using the enclosed AC/USB adapter

Charging the Device

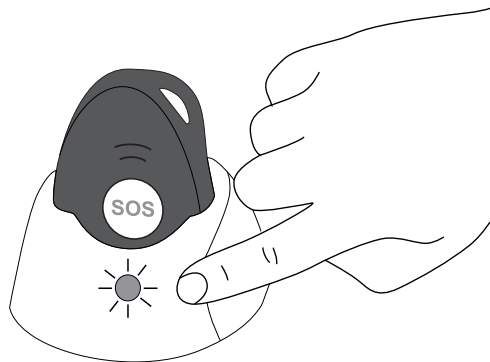

Place the device into the cradle. Rock it slowly until it snaps into place and the red LED light glows.

Turning Device "On"

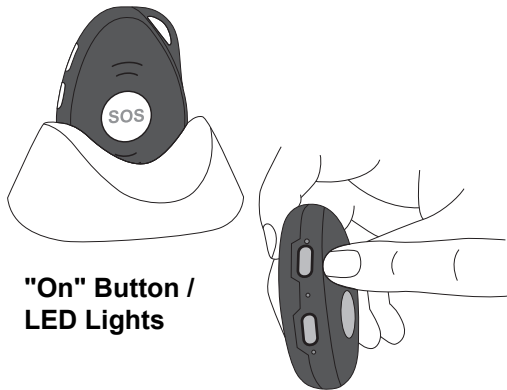

"On" Button / LED Lights


There are three ways to turn the device on:

1. Snap into the charging stand until the red light on the front lights up.
2. Attach the magnetic charging cord.
3. Press the top grey button.

When the device turns on, the LED lights on the side will begin to blink.

Wait at least one minute before pressing the "SOS" button to allow the device time to acquire the cellular network.


Placing a Call

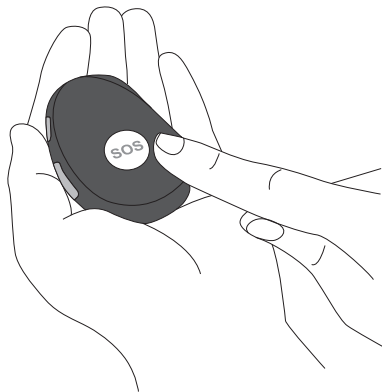

Press the "SOS" button for 2-3 seconds until the device announces "You have activated..." then release. There will be a brief silence while the device places your call

DO NOT PRESS THE SOS BUTTON A SECOND TIME AS IT WILL DISCONNECT THE CALL


Using the Speakerphone

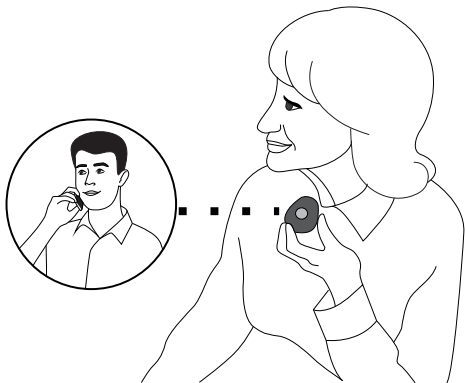

The device will call/text the first pre-set contact. If no answer, it will call the second contact, third, etc. When the call is answered, simply speak to the contact using the built-in hands-free speaker phone.


To adjust the speakerphone volume louder or quieter, press the +/- grey buttons on the side.

Disconnecting a Call

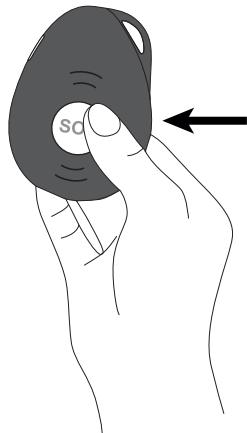


Press the SOS button to disconnect the call

Turning Device "Off"

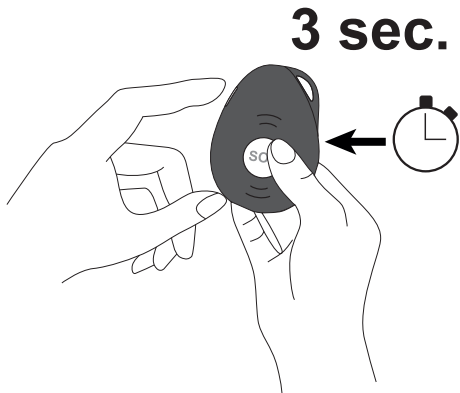

To turn off the device: Press and hold the top side button and SOS button **together** until the device vibrates then let go. The device will power down and the LED lights will turn off within 60 seconds.

Turn Voice Prompts On/Off

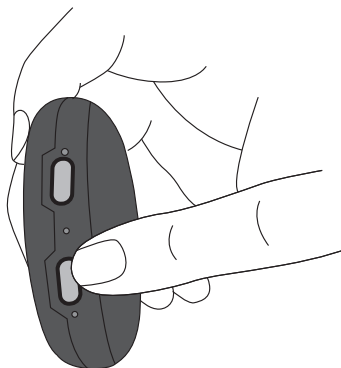

To turn voice prompts On/Off: Double click on the lower side grey button.