



Full mouth, monolithic zirconia reconstruction by Akito Kani, RDT. (Organ Dental Technology)



Recommended for use with Multi-5 **TANAKA ZirColor™ Shift**

Easily and quickly shift your Multi-5 crowns to any shade by dipping them in TANAKA ZirColor Shift Dipping Liquid, designed specifically for Multi-5 zirconia discs.



Multi-5 was designed by Dr. Asami Tanaka, CDT. NMD., world's leading ceramist and a pioneer in zirconia restoration. He is also the founder and president of Tanaka Dental Enterprises.



TANAKA DENTAL USA
8170 McCormick Blvd,
Wi-Fi Professional Building 111
Skokie, 60076 IL.USA
Toll Free: 1-800-325-5266
TEL: 1-847-679-1610
FAX: 1-847-674-5761
info@tanakadentalusa.com
www.tanakadentalusa.com

TANAKA DENTAL Europe
Max-Plank-Strasse 9
D-61381 Friedrechs Dorf
TEL: +49(0)61 72- 8 30 26
FAX: +49(0)6172-84179
service@tanaka.de
www.tanaka.de

ATD Japan
Sugamo 1-14-13
Toshima-ku, Tokyo
Sumishichi Bldg 3F
TEL: 03-5319-1581
FAX: 03-5319-1582
info@atdjapan.co.jp
www.atdjapan.co.jp

TANAKA EnamelZR™


Multi5

Multilayer Zirconia Disk



The Anatomy of Multi-5

Multi-5 zirconia discs feature 5-layers of pre-shading with 3 levels of translucency to produce a smooth, natural gradation.

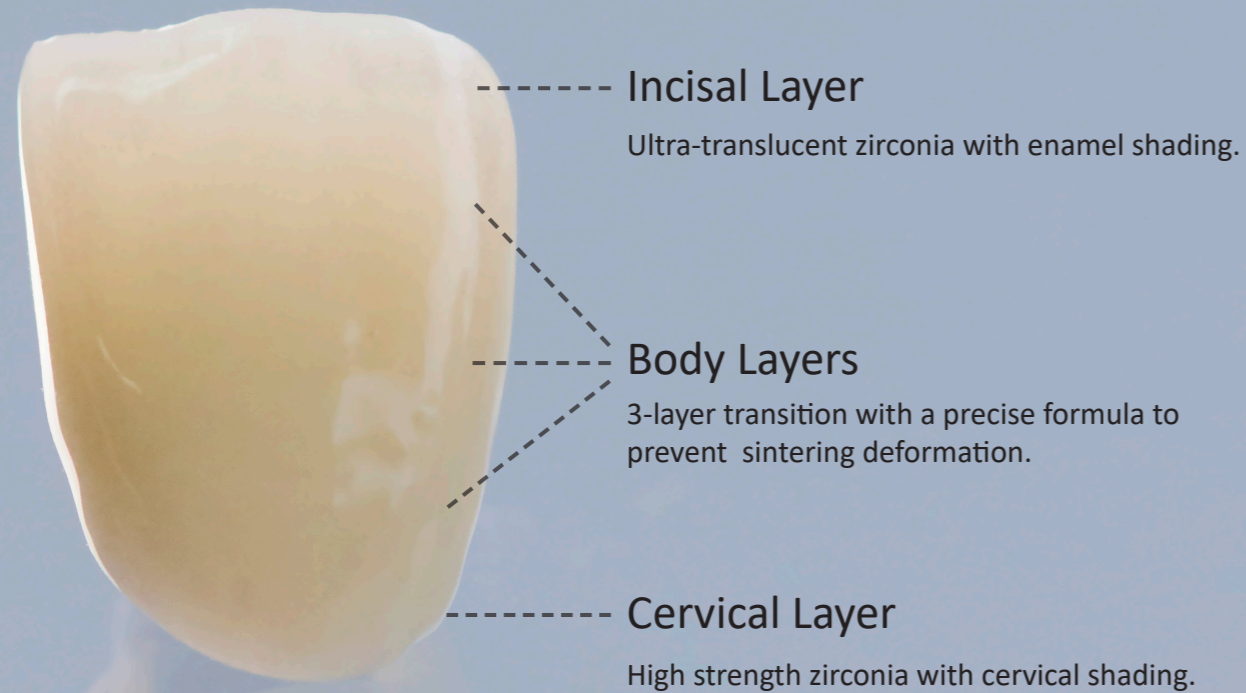


The Difference is in the Enamel

The incisal layer of Multi-5 discs feature a 49% transmittance rating* and a specially formulated enamel shade, which when combined with the other translucent layers, results in a truly natural enamel esthetic.

49%
TRANSLUCENCE

45%
TRANSLUCENCE

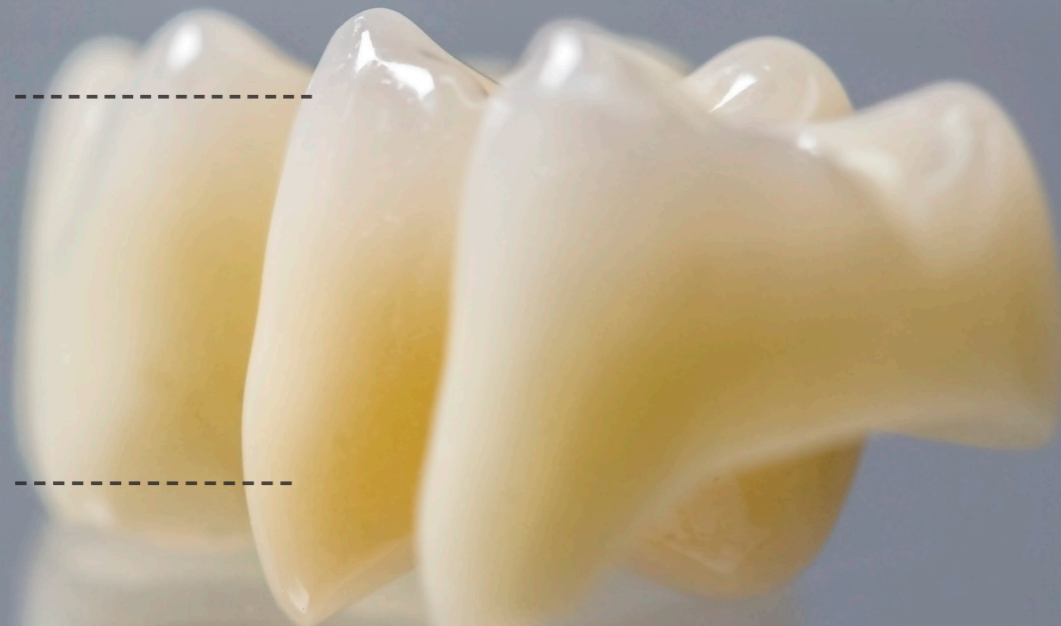

*Evaluated by raw material color.
Wave length of light : 700nm Thickness of sample : 1mm

Restorations fabricated by Yuki Chino, RDT.

Translucent, yet Strong

Our unique layering formula combines high strength layers with super translucent layers to produce zirconia discs that are as strong as they are beautiful.



1100
MPa

12-unit full zirconia, monolithic restoration using Multi-5 by Yuki Chino, RDT.

Multi-5's Wide Range of Application

- ✓ Laminate veneer
- ✓ Inlay/Onlay
- ✓ Anterior Crown
- ✓ Posterior Crown
- ✓ Posterior 3-unit
- ✓ Anterior 3-unit
- ✓ Longspan
- ✓ Full-arch



TANAKA EnamelZR™
Multi-5: Multilayer Zirconia Disk
Flexural Strength: 1100MPa
Translucency: 45~49% transmittance at 1mm thickness
Shade: A1 / A2 / A3 / A3.5
Size: 98x14mm / 98x16mm / 98x18mm
100% Tosoh powder