

User Information Sheet

RONCO Item No : 6335

PRO-TEC **Brand Name**

: Medical N95 Respirator, Vertical Folded Description

Made in Canada

: Non-woven Polypropylene and polyethylene Material

Colour

Characteristics : • Particulate Filter Efficiency in accordance with Health Canada and NIOSH 42 CFR Part 84 (TEB-APR-STP-0059) evaluation: ≥95%

• Breathing resistance/Differential Pressure: a. Inhalation: <35 mmH₂O (343.2 Pa); b. Exhalation: <25 mmH₂O (245.2 Pa)

Bacterial Filtration Efficiency (ASTM F2101-19): 100% BFE


Ability to withstand 10 N (0.98 kgf) at connection points of each headstrap

• Pass Fluid resistance (ASTM F1862/F1862M-17): Level-3 Pass at 160 mmHg (21.3kPa)

• Biocompatibility (each skin contacting component): Pass the Cytotoxicity, Skin Sensitization and Skin Irritation tests

Standards : Meets Health Canada guidelines for N95 respirators.

Specification

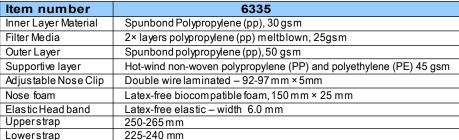

: 30 masks x 6 boxes = 180 masks per case Packaging

Product Features

- : The lightweight, vertical folded design of the respirator provides reliable, and convenient worker respiratory protection especially in the healthcare environment
 - Individual automated packaging helps protect the respirator from contamination before use
 - · Compatible with a variety of eyewear
 - Conforms to eye and nose contours of the face ensuring more room for eyewear
 - Adjustable noseclip with soft and comfortable fog-reducing foam
 - This respirator contains no components made from natural rubber latex
 - High-quality noseclip ensures proper fit to all facial contours
 - Headstraps are 100% latex-free to minimize risks of allergic reactions
 - Excellent resistance against bodily fluids and blood
 - One Size Fits All

Applications : • Operating Rooms

Clinics

Laboratory

TB Wards

Emergency rooms

· Pandemic preparedness planning, stockpiling

· Patient Care

• Grinding, Sanding, Sweeping, Bagging and other dusty or arid operations

Labor and Delivery

• Can also be used to help reduce inhalation of certain airborne biological

 Infection Control Practices particles like mold, Bacillus anthracis, Mycobacterium tuberculosis

Storage & Handling

: Store in original packaging in a cool,dry and well-ventilated area, away from dust, direct sunlight, moisture and excessive heat (-20 to +30°C) (-4 to +86°F). Avoid direct contact with the used mask, Put it in a plastic bag and discard used masks into a proper lidded trash bin.

Instructions & Restrictions For Use:

- Before use, the wearer must read and understand the User Instructions provided as a part of the product packaging.
- This product is for single use and should be discarded after each use or a maximum of 8 hours of continuous use.
- Do not use in atmospheres containing less than 19.5% oxygen.
- Do not use in atmospheres immediately dangerous to life or health (IDLH).
- Do not exceed maximum use concentrations established by regulatory standards.
- Failure to properly use and maintain this product could result in injury or death.
- All approved respirators shall be selected, fitted, used and maintained in accordance with OSHA, CSA and other applicable regulations. This respirator does not have any replacement parts.
- This respirator helps protect against certain particulate contaminants but does not eliminate exposure to or the risk of contracting any disease or infection.
- Misuse may result in sickness or death.





