



THE LEADER IN VISION DIAGNOSTICS™

E-Z TILT CHAIR



Attractive, modern design.



Sleek, low profile hydraulic base is more attractive and less obtrusive.



Misty-gray painted base blends naturally with overall color scheme.



Fast reclining motion. The patients' own weight counterbalances the reclining motion throughout the full range of travel.



Exclusive double hydraulic cylinder system provides an extended 300mm range of motion — from extra low to extra high.



Convenient single-handed and foot operated controls for chair reclining, elevation, rotation, and headrest adjustments.



Unique headrest locking indicator insures patient safety and stability.



It Doesn't Get Any Easier.

The **MARCO E-Z TILT CHAIR** combines the most simplified manual reclining motion with a new modern, sophisticated appearance. Utilizing an effortless, silky smooth counterbalancing system, the E-Z Tilt makes reclining any size patient easy and instantaneous; no electronic motors to slow you down! The Marco E-Z Tilt also incorporates the only double hydraulic cylinder system in the industry, creating an extended range of elevation — from extra low to extra high, that accommodates a full range of patient and practitioner sizes. The sealed hydraulic chair base is an attractive, low-profile non-leaking system that automatically “bleeds” unwanted air bubbles from the hydraulic hose. This exclusive design also



Exclusive telescoping hydraulic cylinder system.



Ergonomically designed arms rotate up for patient flexibility.



Foot-activated rotation lock can easily be activated from both sides of the chair.



Single-handed headrest with locking indicator.

eliminates many complicated and unnecessary components, resulting in much greater reliability.



New overlay switch panels on both sides of the chair for convenient fingertip adjustments.



MEASUREMENTS

Ask your Marco Distributor about The Marco Package Plus Program.



PLAY IT SAFE.

Protect your warranty. Buy your Marco product only from an authorized Marco distributor. Call Marco for the name of the distributor nearest you.

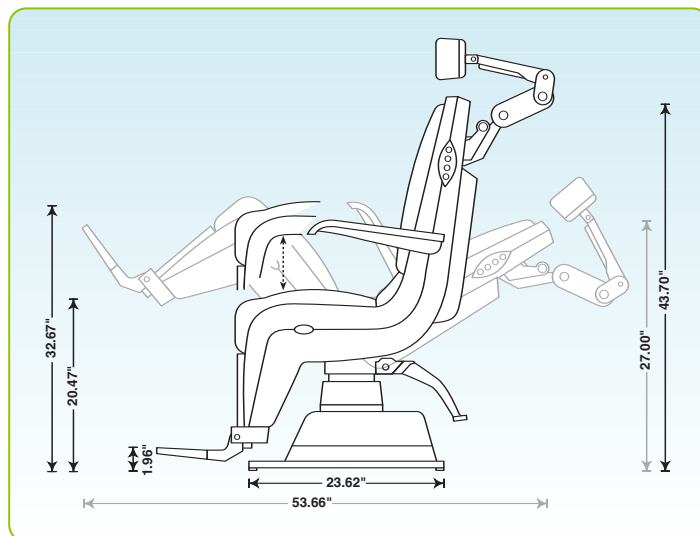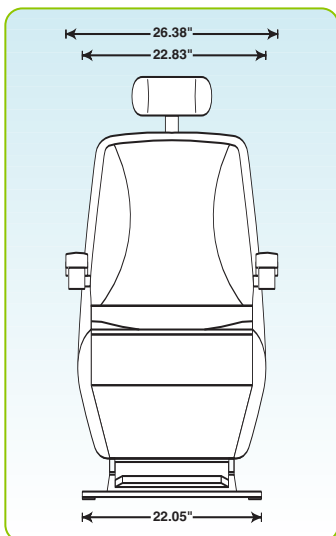

PO55000 BR-EZ

T H E L E A D E R I N V I S I O N D I A G N O S T I C S™



11825 Central Parkway | Jacksonville, FL 32224 | 800 874-5274 | 904 642-9330 | www.marco.com