

Kowa



Technology for Life Science

KOWA APPLANATION
TONOMETER

HA-2




Convenient Item For Daily Clinical Practice



Technology for Life Science

KOWA APPLANATION TONOMETER HA-2

● Standard set of HA-2

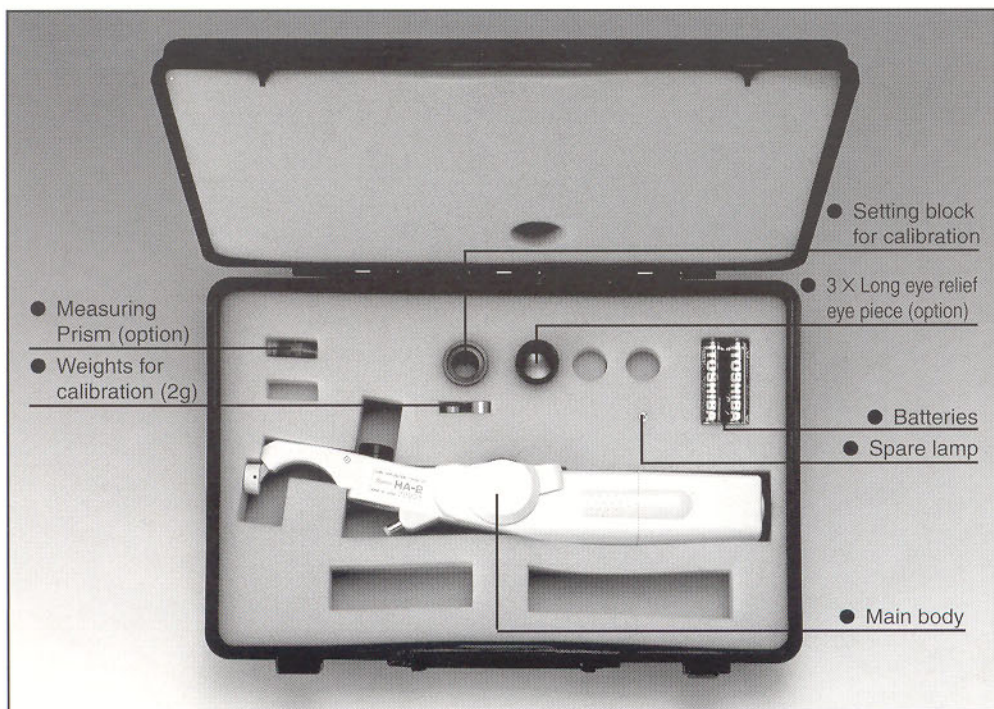

3 X Long eye relief eye piece (option)

● Main Features

- ◆ Compact and light-weight, well balanced and easy to handle by single hand.
- ◆ Unique one-spring mechanism assures correct reading regardless of patient posture.
- ◆ Built-in illumination permits you to use the instrument regardless of the illuminating condition of consultation room.
- ◆ Direct reading type precision dial requires no conversion table to change the readings.
- ◆ Accessories such as adjustable head rest and interchangeable eyepiece provide most suitable conditions for measurement.

● Specifications

Type	Hand-held applanation tonometer
Applanation prism	Doubling prism interchangeable (option)
Scale range	0~60 mm Hg (1 mm Hg division)
Illumination	
a) ON/OFF	Interlocked with thumb wheel
b) Bulb	Lens bulb 2.5 V, 0.25 A interchangeable
c) Battery used	1.5 V, 2 pcs Interchangeable
Eyepiece	
Magnification	6 X
Head rest	Movable
Length	290 mm
Weight	240g with batteries

* All specifications and external design are subject to modification of change without prior notice.



Kowa Kowa Company, Ltd.

World Sales Headquarters
 4-14, Nihonbashi-honcho 3-chome, Chuo-ku, Tokyo 103-8433 Japan
 Phone: 81-3-3279-7331 Facsimile: 81-3-5255-7516
 URL: <http://www.kowa.co.jp/e-life/>

Hamamatsu Factory
 3-1, Shinmiyakoda 1-chome, Hamamatsu City,
 Shizuoka Pref., 431-2103 Japan

Kowa Europe GmbH Kowa Optimed, Inc.

Immermannstrasse 43B
 40210 Düsseldorf F.R.Germany
 Phone : 49-211-1793540
 Facsimile : 49-211-161952
 URL : <http://kowa-europe.com/>

20001 S. Vermont Ave. Torrance,
 CA 90502, U.S.A.
 Phone : 1-310-327-1913
 Facsimile : 1-310-327-4177
 URL : <http://medical.kowa.com/>



This catalog was made with 100% recycled paper

Printed in Japan 0594500MA-2