

NG Feeding Tubes (Yellow)



VS

NG Plus Feeding Tubes (Purple)

Tip is two sizes larger than the rest of the tube and only has 1 side hole



Moderate kink resistance

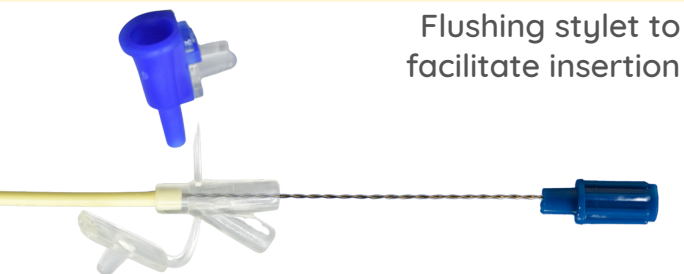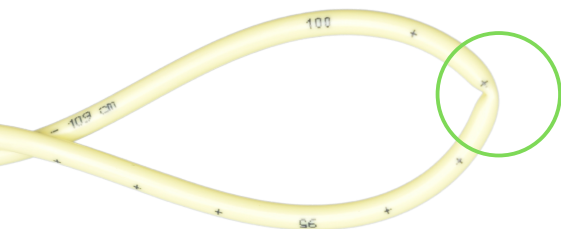

Flushing stylet to facilitate insertion



Y-adapter for feeding, medication and flushing


Yellow, opaque, radiopaque polyurethane with depth markings




Uniform diameter with 4 side holes




Superior kink resistance when exiting the nostril

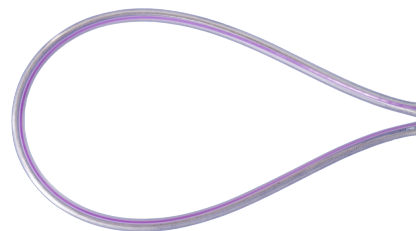

Flushing stylet to facilitate insertion AND designed to allow easy wire removal in tight turns

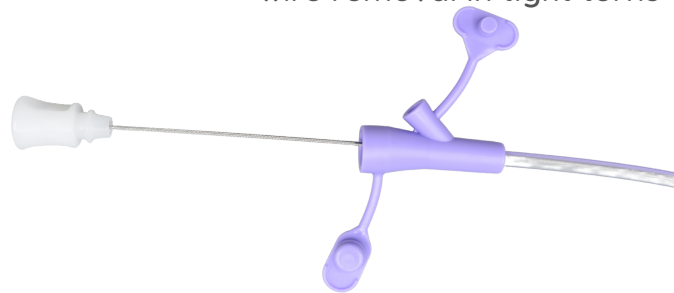

Dual port designed for funnel syringe and luer slip syringe or connect directly to feeding tube extension set

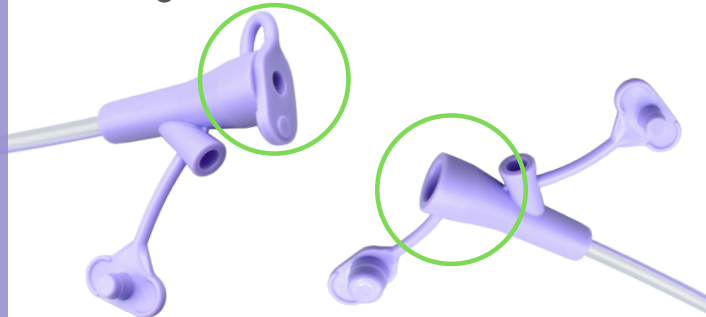

Clear, LOW FRICTION polyurethane with radiopaque purple stripe and depth markings



Color coded purple for safety to indicate feeding and prevent accidental misconnections