


SHORT HANDLE EXTRACTION KIT





This kit consisting of a short handle with different luxators is very useful for the extraction of wolf teeth, incisors, canines and deeply located root fragments. It contains straight, concave and convex curved luxators in a 4 and 5 mm version. Particularly the concave curved tips allow a very good access to the palatal side of the wolf tooth and enable a careful extraction with minimal damage. The short handle is designed to allow the gentle usage of a hammer as well. Furthermore, this short handle can be combined with all different inserts from the Dental Pick System.

Shapes

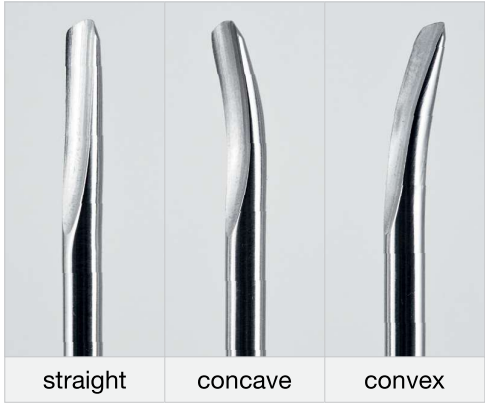

Widths

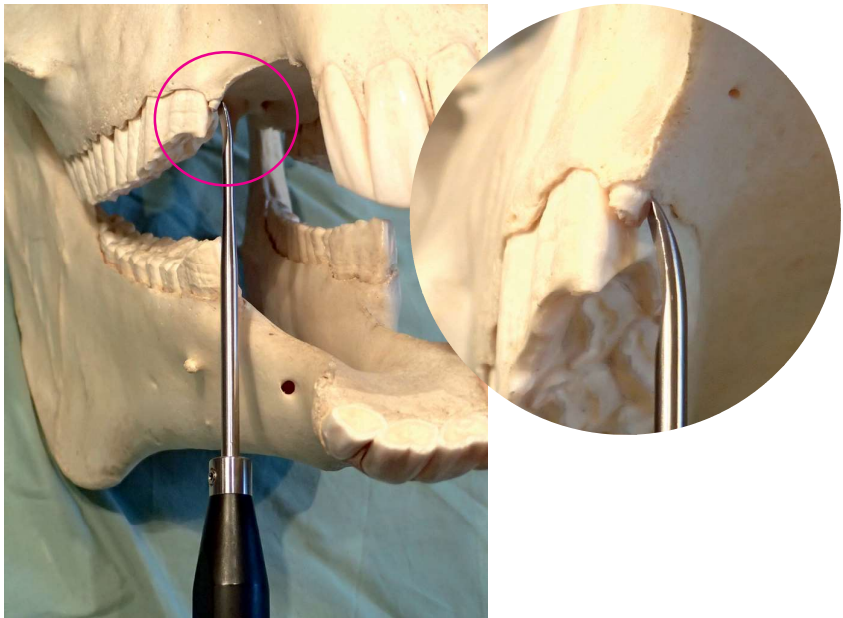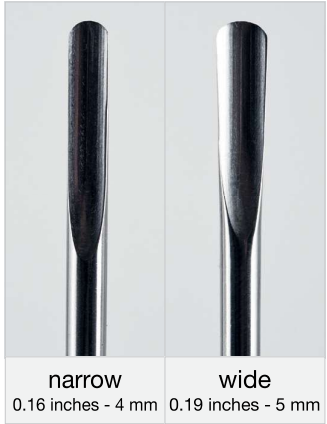

Art.-No. #	Description
PEGA_312001	Handle (short, straight)
PEGA_312004	Luxator Insert 150/4 straight
PEGA_312002	Luxator Insert 150/5 straight
PEGA_312005	Luxator Insert 150/4 concave
PEGA_312003	Luxator Insert 150/5 concave
PEGA_312006	Luxator Insert 150/4 convex
PEGA_312007	Luxator Insert 150/5 convex
PEGA_312040	Wolf Tooth Extraction Kit Contains a short handle and 4 luxators (straight and concave curved inserts in 4 and 5 mm width)
PEGA_312041	Extraction Kit with short Handle and 6 Luxators Contains a short handle and 6 luxators (straight, concave and convex curved inserts in 4 and 5 mm width)