


ISSUE DATE:	18-JAN-2021	REVISION:	001	DOCUMENT NO.	MA-SPEC-3P
-------------	-------------	-----------	-----	--------------	------------

AMAX LEVEL-3 (PINK) PROCEDURAL FACE MASK SPECIFICATION

Product Code: AMSCL3P

LOT#: _____

Description: 6.9" x 3.75" (Level 3 Mask Pink with Aluminum Wire and Ear-Loops)

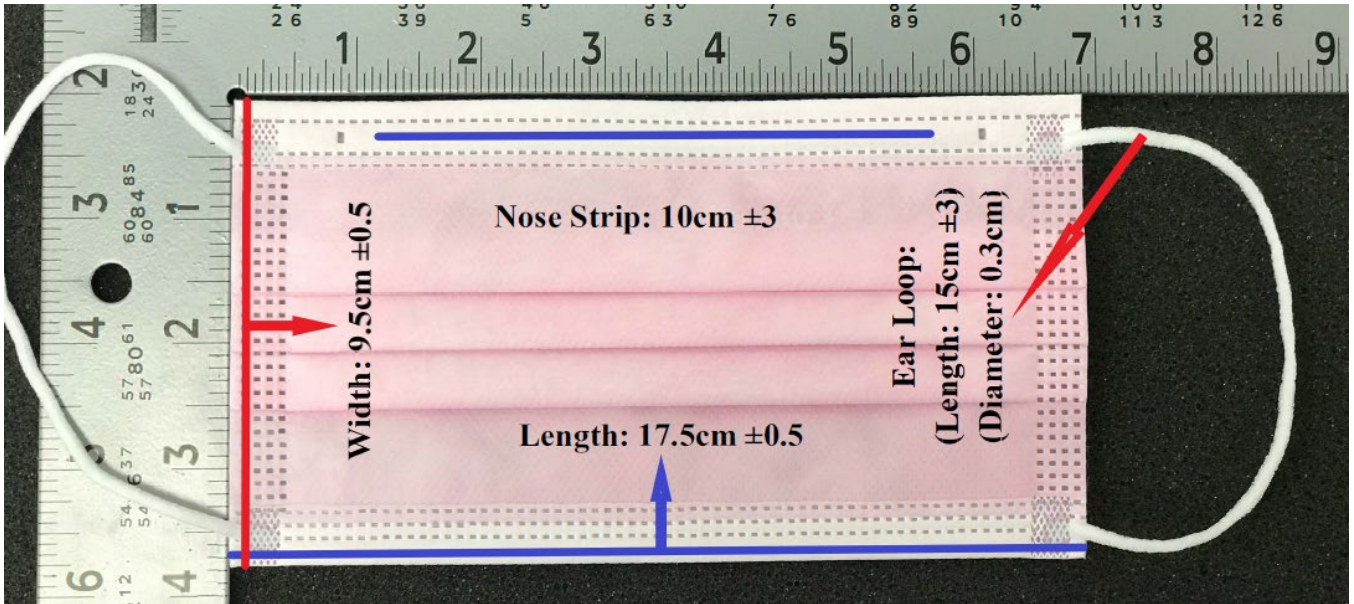

GS Medical Confidential

Item No.	Body Size		Inner Layer Spun-bond non-woven fabric Grade SSS White 25 gsm	Middle Layer Melt-blown non-woven fabric High Filtration White 20 gsm	Outer Layer Spun-bond non-woven fabric Grade SSS Pink 25 gsm	Nose Strip Aluminum Wire 10cm ±3	Ear Loop Elastic Polyester and Spandex Length: 15cm ±3 Diameter: 0.3cm
	Length	Width					
AMSCL3P	17.5cm ±0.5	9.5cm ±0.5					

QC Use Only

Average face mask ear loop tensile strength must be 10N minimum.

Label Issuance and Reconciliation				
Label Type	Quantity	Date of Issue	Quantity Remaining	Disposal

Production Summary

Production Start Date: _____

Production Finish Date: _____

Total Pallets: _____

Total Cases: _____

Total Units: _____

Operator and Packers Use Only

Line set-up and clearance sign off: _____ Date: _____

Packer sign off on line clearance: _____ Date: _____

Information check for Inner and Outer labels (Packer Sign off)

I have read the label and all information pertaining to this product is correct.

Signature: _____ Date: _____

Inspection Components:

- Check Mask quality of 3 random mask samples every 15 minutes.
- Visually inspect fabrics ultrasonic welds with no loose end.
- Ensure mask ends are neatly trimmed and Ear Loops are firmly attached.

Approved by (QC): _____ Date: _____