

# Foot Spa

## Operating Instructions



# Operating Instructions

## Specification

Model: Suk-803  
Input: Ac 110V

Power: 60W  
Output: 0-15V; 0-2.5A

Freq: 50Hz-60Hz

## Getting Started



1. Before using, plug the power cord, array and wrist strap into their responding jacks, and then pour some tap water into the foot basin (please note that this step can't use distilled water, or the system won't work). Afterwards, you need to wear the wrist band, submerge the array as well as your feet in water and then turn on the power. **Make sure the level of water goes over the top of the array.**
2. **Functional Screen:** The unit is functioning when the display screen flashes as follows after turning on the power switch:

```
Enter Time 10—30  
Time=30 V3.0
```

- **Select operating time:** By default, the Foot Spa is under the 30-minute working cycle. To adjust this, use the + and - or buttons to increase and decrease the value.

```
Time=06: 28 pos+  
T=30 M=3 E=1.0
```

- **Select operating mode:** By default, the Foot Spa is under Mode 3. Select a different mode by pressing + or - buttons. Once the working period is done, the machine will stop automatically.



3. It may be necessary to add a pinch of salt to the bath water to meet the desired current output and enhance ionic conduction. (**Mode1: 0.6 amps, Mode2: 0.8 amps, Mode3: 1.0 amps. If the concentration is too high, the machine will beep. You may pour some tap water or select a low mode**)
4. This product is the third generation product and equipped with a function of recording frequency as well as operating time. Examine the record by simultaneously pressing 1 and ENT. For example:

00298	Mins
93	Time

The first line means the total operating time and the second line refers to the frequency. When the operating time reaches to 10 minutes, the frequency will increase 1 time. Back to the menu by pressing ESC.

Wipe the data by simultaneously pressing 2 and ENT. The screen reads:

PAS: V2.0
Clear Contents

Press the default password 666666 and then presses (ENT) the wipe the records. Afterwards, the screen reads:

Process: OK
Mins=00000

Press ESC to go back to the menu.

5. Change the default password by simultaneously pressing 3 and ENT.

The screen reads:

PAS: V4.0  
Change Password

At first, input the default password and ENT, the screen reveals:

Process: OK  
New PAS: 123456

Second, input a new six-digit password and press ESC.


6. The Far Infrared Charcoal Bamboo Belt can be operated in conjunction with the Foot Spa, or by itself. Place the belt on your body. Most people wrap it around their lower back. You may wrap your leg with the elastic band. Turn on the corresponding power switch, select desired mode-H, M or H and the belt is now working. The buttons that read L, M and H correspond to Low, Medium and High. Be sure to hit the power switch again at the end of the session to turn the belt off.



**Control panel for the far infrared charcoal bamboo belt.**

## Attention

1. Do not connect or disconnect the plug to a household outlet with a wet hand, please pay attention to your personal safety.
2. Avoid striking the machine.
3. Wipe this product with a cleaning cloth. Do not clean the machine with solvents, such as benzine, thinner, or commercially available cleaners.
4. Gently press the buttons to ensure continued use.
5. Contact the wrist strap before use-make sure it is secure, but firmly attached.
6. You may lift up the main body with the handle while the machine is working for your convenience.(Spinning the handle after pull it out the main body slightly)

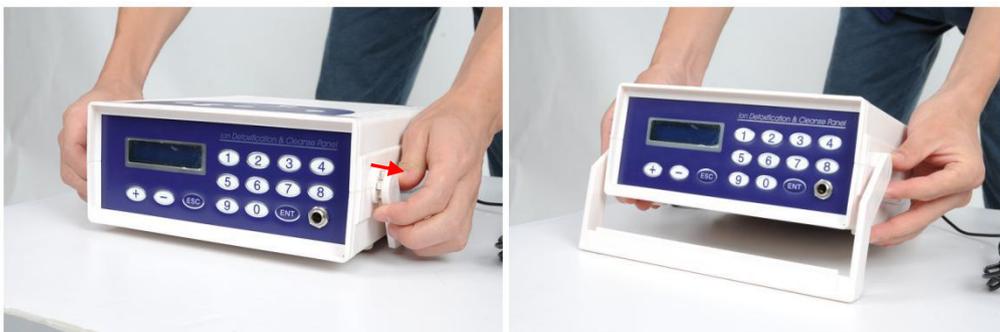

## Troubleshooting

### Is the LCD screen not lighting up?

Make sure that this unit is powered on. What's more, make sure there is no fuse blown.

### Is the array stop working?

Make sure that cable is securely connected. Wipe the array with a cleaning cloth from the center outward to the edge. If the above two steps fail, you need to replace the array. (The array may be used for 50 to 60 sessions.)

### Unresponsive after pressing a certain button?

Turn off power switch and then on in about 10 seconds.

## Packaging List

Carrying Case	Main Body	Array	Power Cable	Wrist Strap	USB Cable	Far Infrared Belt	Manual
1	1	1	1	1	1	1	1

## Notice

1. This machine must **not** be used when you are:
  - ◆ A pregnant woman or breastfeeding mother;
  - ◆ With a type of implanted electromagnetic device;
  - ◆ An organ transplant recipient;
  - ◆ A child younger than 8 years old.
2. People who have hypoglycemia need to feed themselves before using this unit.
3. Please do not immerse the wound under water.

