

Mistral-Air®

Forced Air Warming



Fast
Safe
Easy

Warm Air with Passion

Mistral-Air® warming blankets are very soft, strong and convenient. The diffusing surface ensures optimum performance and increased heating capacity. It also enables the use of low pressure, preventing the blanket from floating and allowing it to drape perfectly.



Mistral-Air® Forced Air Warming System



**A Fast,
Safe & Easy
System**

Blower / Blanket tuned design

The Mistral-Air® warming blankets and Mistral-Air® Plus warming unit are attuned to each other. This ensures optimum performance and equal distribution of heat throughout the entire blanket



Mistral-Air® Plus Warming Unit

The Mistral-Air® Plus warming unit is very powerful, silent and provides a constant temperature to your patient within only 60 seconds. The default start up temperature of 38°C ensures the immediate start of therapy and reduces the risk of human error.

Safety first

In addition to the over-temperature alarm, the under-temperature alarm is a standard built-in safety factor guaranteeing that the temperature which has been set is delivered to your patient.



Mistral-Air® Blankets

The wide range of Mistral-Air® warming blankets ensures there is a blanket for each application and can be used for pre-, per- and post-operative warming.



Mistral-Air® Reflective Blankets

This range of warming blankets is equipped with a reflective top layer, increasing heating capacity and reflecting your patient's radiant body heat, preserving body temperature.



Mistral-Air® Koala Quality Blankets

This range of warming blankets consists of under bodies which allow you to place your patient on top of the under body. These warming blankets are equipped with a reflective layer and have adhesive strips attached so you can "tuck in" your patient. Like a young Koala kept warm and protected in its mother's pouch.

More products from TSCI:
Fluido®*AirGuardSystem*
High Flow Blood & Fluid Warming



Fluido®*Irrigation*
Fluid Warming



Thermoflect®
Heat Reflective Technology



Partner in Patient Temperature Management

TSCI is a leading European company in the field of hypothermia and offers a complete range of innovative solutions for patient temperature management.

Fluido® and Mistral-Air® are registered trademarks of TSCI.

Fluido® products are covered by the following European patents:1313521 and 1446179.

Other patents pending.

Thermoflect® is a registered trademark of Encompass Group, LLC.

TSCI

The Surgical Company
International B.V.

Beeldschemweg 6F
3821 AH Amersfoort
The Netherlands

T +31 (0)33 450 72 50
F +31 (0)33 450 72 60
contact@tsci.nl
www.tsci.nl

 The Surgical
Company
GROUP


MA0220 / 320
Adult




MA0230 / 330
Paediatric




MA0250 / 350
Lower Body

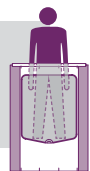


MA0260 / 360
Upper Body

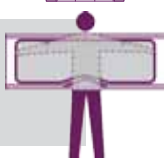

MA0265 / 365
Half Upper Body

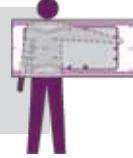

MA0270
Torso

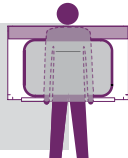


MA0280 / 290
Surgical Access



MA0510
Tube



MA0340
Neonatal

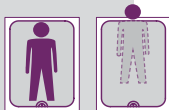

MA0400
Full Underbody



MA0450
Underbody



Mistral-Air® Forced Air Warming Blankets

Mistral-Air® Blankets

Article number	Description	Dimensions
MA0220	Adult	2180 x 1270mm
MA0230	Paediatric	1450 x 1020mm
MA0250	Lower Body	1420 x 1020mm
MA0260	Upper Body	2000 x 760mm
MA0265	Half Upper Body	1600 x 760mm
MA0270	Torso	760 x 1240mm
MA0280	Surgical Access (sterile)	2200 x 1270mm
MA0290	Surgical Access	2200 x 1270mm
MA0510	Tube	1750 x 150mm

Mistral-Air® Reflective Blankets

Article number	Description	Dimensions
MA0320	Adult	2180 x 1270mm
MA0330	Paediatric	1450 x 1020mm
MA0350	Lower Body	1420 x 1020mm
MA0360	Upper Body	2000 x 760mm
MA0365	Half Upper Body	1600 x 760mm

Mistral-Air® Koala Quality Blankets

Article number	Description	Dimensions
MA0340	Neonatal	850 x 760mm
MA0400	Full Underbody	2150 x 1020mm
MA0450	Underbody	1630 x 1020mm

TSCI

The Surgical Company
International B.V.

Beeldschermweg 6F, 3821 AH Amersfoort, The Netherlands. T +31 (0)33 450 72 50, F +31 (0)33 450 72 60

contact@tsci.nl, www.tsci.nl VAT No. NL 8111.49.948.B.03 Chamber of Commerce (Amersfoort) No. 31029240

Mistral-Air® Forced Air Warming



Mistral-Air® Plus Warming Unit

Article number	MA1100-EU
Voltage	220-240VAC
Frequency	50/60Hz
Current	3A
Average power	525W
Dimensions	276mm x 385mm x 239mm
Weight	± 6.0kg
Hose length	1.8m
Power cord length	4.0m
Set temperature	32°C, 38°C, 43°C or Ambient Air
Filtration	HEPA, 0.3µm, 99.99%, H13 conform EN 1822
Limited Warranty	3 years
Classification 93/42/EEC & 2007/47/EC	Class IIb
Classification IEC 60601-1	Class I, Body Floating (BF)
Classification IEC 60529	IP21



Mistral-Air® Pole with Basket

Article number	MA5100-B
Dimensions pole	1000mm height, ϕ 25mm
Dimensions Base	ϕ 600mm
Castors	5, anti-static, latex free, 2 brakes
Basket dimensions	230mm x 230mm x 250mm
Weight	± 6.7kg



medical
sales &
service

A: 38 Hector Street OSBORNE PARK WA 6017
T: 08 9361 0843 | **F:** 08 9470 3647
E: medsales@medsales.com.au | **W:** www.medsales.com.au