

seca 762

Mechanical personal scales
with 500 g graduation



- Robust construction.
- Large, easy-to-read dial.
- Low profile and wide platform.
- Precise and durable.

seca 762:

Classic highlight in weighing technology.

This seca classic has greatly influenced the image of scales. With an absolutely precise dial, synthetic covering used in medical applications, and steel plate housing from the good old days. The modern age has also left its mark on these scales, which are available in many different colors and models. For one thing a classic should never be boring.

Modern classics for several decades.

The precise design classic has a large clear and well defined dial. With a large, low-level, non-slip platform and scratchresistant finish it is easy for the patient to step on and weigh.

The face of the dial is especially large to make reading the figures much easier.




The rugged metal housing withstands years of use under the hardest conditions imaginable.

seca 762

Technical Data

- Capacity:
150 kg, 320 lbs, 150 kg / 320 lbs
- Graduation: 500 g, 1 lbs, 500 g / 1 lbs
- Dimensions (WxHxD):
305 x 120 x 470 mm / 12 x 4.7 x 18.5"
- Weight: 3.5 kg / 7.7 lbs
- Functions: reset-to-zero function



seca clara 803

The rediscovery
of clarity.

product sheet { **seca clara 803**



- Two-component rubber coating. For pleasant footing.
- Four stable rubber feet. For a safe standing.
- 28-millimetre LCD figures. For easy reading.
- More than 165 years of experience. For maximum precision.
- From the global market leader for medical weighing and measuring.

seca clara 803:

The rediscovery of substance.

First take a look at the clear design and then at the finishing: even our simplest scale is complex in its design and robust in its construction. It gives even people who weigh up to 150 kilograms safe footing.

The rediscovery of feeling good.

Go ahead, pick up our starter scale and take a close look. You'll be impressed by its classic design and its surface gleaming with high-quality materials. The high-quality rubber coating consists of a special combination of two components which lends softness and solidity. The material can stand up to all commercially available cleaning and disinfecting products.

Professional can be so private.

Whether black or white, whether for professional or private use, this seca starter scale satisfies the highest demands made by private persons and anyone who appreciates high-quality workmanship. The seca clara 803 has its strengths on the outside and on the inside. Underneath the weighing platform you'll see longitudinal and cross braces set one centimetre apart, giving the scale unparalleled stability.



Lightening the load.

Thanks to its low tare weight of only 1.4 kilograms and the handle built onto the back, the scale is easy to carry, making it suitable for users on the move. Viewing the measurements is just as easy from the unconventionally big display with 28-millimetre-high figures. You see, seca cuts no corners even with its lowest-priced scale.

28-millimetre LCD figures. For easy reading.



Close ribbing for a sure foothold.



seca clara 803

Technical data

- Capacity: 150 kg / 330 lbs / 23 sts
- Division: 100 g / 0.2 lbs
- Dimensions scale (WxHxD): 315 x 35 x 325 mm
- Weight: 1.4 kg
- Functions: Auto-HOLD, automatic switch-off
- Parts to be delivered: 4 AAA batteries, 4 extension feet for carpets

seca robusta 813

The standing,
taking everything lightly.



- Solid steel frame. For persons up to 440 lbs.
- Non-slip rubber base. For safe standing during measuring.
- Big LCD figures. For easy reading.
- Step-off function. For weighing without switching on scales.
- More than 165 years of experience. For maximum precision.
- From the global market leader for medical weighing and measuring.

seca robusta 813:

The standing, making everything look a breeze.

Take a look at the extra-wide weighing platform. Even very big people can stand here safely and securely. People of larger girth too can be accommodated. Our intention is to lighten the task for weighty people.



A scale of standing, making the sturdy look beautiful.

More and more people weigh more and more. Now there's a stable scale for them. Thanks to its solid steel frame, it can withstand 440 lbs and makes light of heavy use over many years. It's also no burden to step onto the scale, since its construction is extremely flat. That's a big advantage for people whose size makes it difficult for them to get about.

The standing, making the robust look elegant.

All rubber is not created equal. For this weighing platform seca has chosen a high-quality two-component nubby material which is soft like rubber and pleasant to the touch. But don't be fooled. This material is also tough. It fulfils two requirements at the same time: pleasant weighing and an easy-to-clean surface that stands up to all commercially available cleaning and disinfecting products.

The standing, paying attention to details.

Weight and girth are no problem on the seca robusta 813. You don't even have to bend over to switch it on, thanks to the unique Step-off function. An additional problem with many scales is the display. Not so here, because thanks to the 1.1 inch-high LCD figures, it is extremely easy to read the measurements. You see, a scale from seca makes life easier by paying attention to details.

High-quality, two-component nubby material feels as soft as rubber.



Big LCD figures. For easy reading.

seca robusta 813

Technical data

- Capacity: 440 lbs / 200 kg
- Division: 0.2 lbs / 100 g
- Dimensions scale (WxHxD):
17 x 1.8 x 14.7 inch
- Weight: 5.7 lbs / 2,6 kg
- Functions: Step-off, Auto-HOLD, automatic switch-off.
- Parts to be delivered:
4 AAA batteries,
4 extension feet for carpets

seca 876

Flat scales for
mobile use

new



- Mother-Child function for trouble-free weighing of small children.
- Large leveling base for high stability.
- Easy to transport – scale itself weighs only 3.6 kg.
- Battery power for sites without electrical outlets.

seca 876:

A lightweight for heavyweight use.

In the seca 876 you've found the right solution if you're looking for a lightweight scale (only 3.6 kg) that's stable, compact, easy to clean and equipped with a Mother-Child function.



Light but stable.

The base is made of a mixture of ABS plastics and Plexiglas. This combination guarantees a material with outstanding surface hardness and impact strength. Furthermore, it is easy to clean with any commercially available disinfectant. Additional steel elements in the base make sure that persons up to 250 kilograms can weigh themselves easily.

Light but impact resistant.

ABS plastics and Plexiglas are not only very stable but are also especially light, a quality that is important for those interested in mobile use of the seca 876. Just as important for use away from home is the large leveling base which allows the scale to be positioned without any trouble regardless of the surface underneath it. It's easy to step onto the scale, thanks to the flat construction and slip-proof standing surface.

Light, but with Mother-Child function.

The scale is easy to operate even with the display and operating elements integrated in the base. Just a light tap of the toe on the scale surface turns on the scale. Infants and small children can be easily weighed while being held by the mother. All she has to do is weigh herself separately first. Then the Mother-Child function ascertains the tare of the weight.



Flat construction and slip-proof surface for a comfortable foothold.



Large leveling base on the bottom of the scale for safe placement.



Large display for easy reading.



Tip-on automatic.


seca 876

Technical Data

- Weight capacity: 250 kg
- Scale division: 100g < 150 kg > 200 g
- Scale dimensions (WxHxD): 313 x 58 x 315 mm
- Weight of scale: 3.6 kg
- Power: batteries
- Functions: 2 in 1 (Mother-Child function), Auto-HOLD, Tip-on automatic for switching ON and OFF. kg/lbs/sts switch-over

seca 874

Flat scale for mobile use

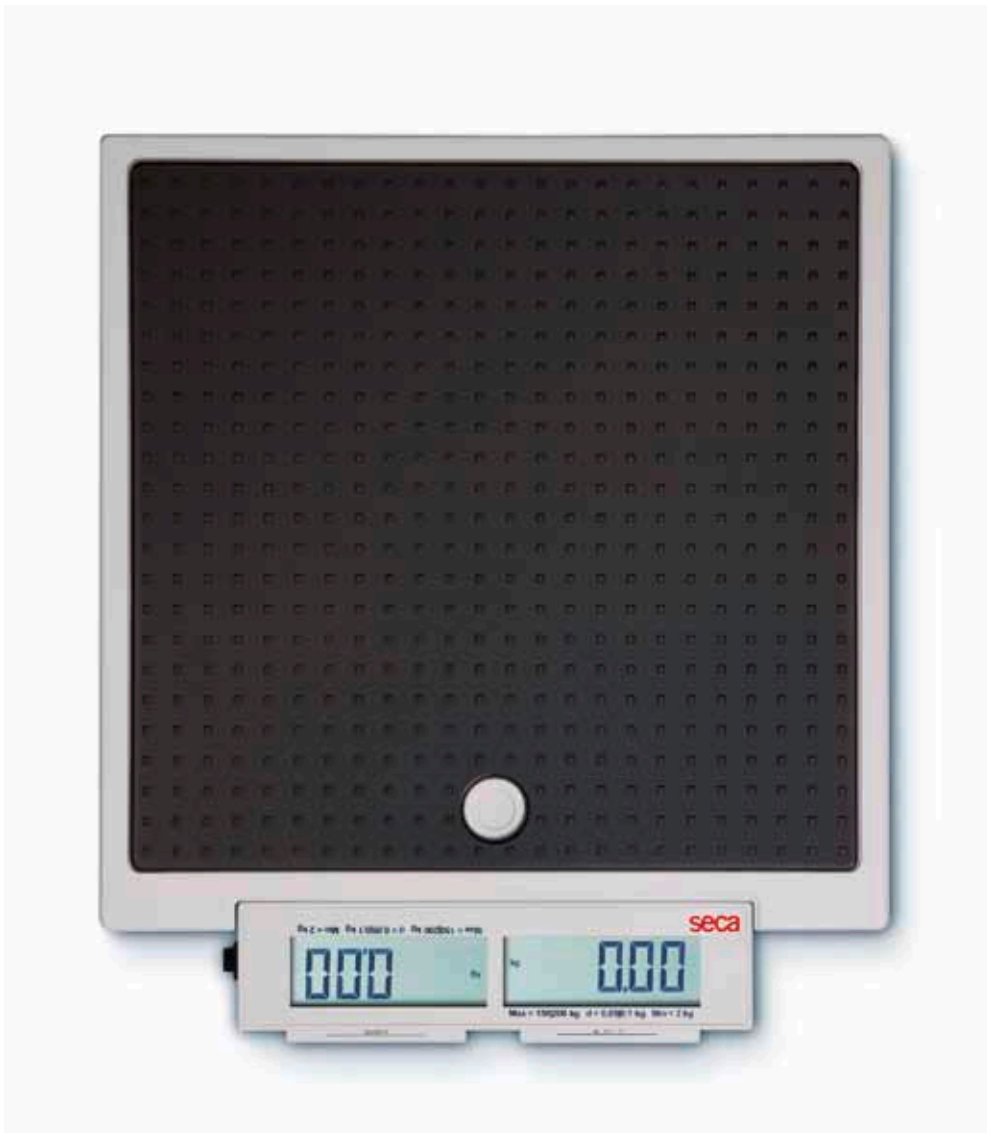

- Mother/child-function for trouble-free weighing of small children.
- Large leveling bases for increased stability.
- Low device weight and battery power for mobile use.
- Double display for readout from different positions.
- Two push buttons for effortless operation with the foot.

seca 874:

The scale that thinks of people.


The seca 874 can do everything that a scale in mobile use should do and then some. It features extremely robust construction and easy-to-use design with long-life batteries and precision unaffected by humidity levels. The scale is so simple that you need only the tip of your toes to operate it.



Great stability.

The seca 874 is very light at only 8.9 lbs, but also very tough. Reinforced by steel elements in the frame, the scale also boasts a top-quality ribbed surface, extra-strength two-component ABS plastic and Plexiglas. Small children and adults up to 440 lbs can be weighed accurately. Especially large leveling bases keep the scale stable on any type of flooring.

Both doctor and patient can see the result at the same time in the double display.

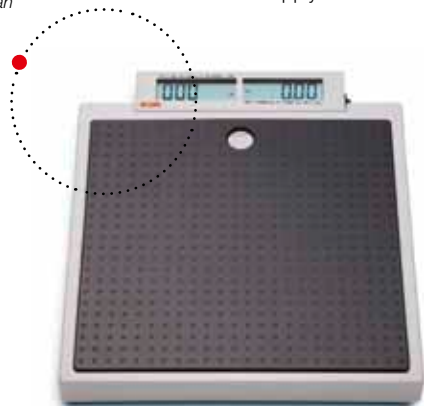

Simple operation.

When conducting one examination after the other, there's no need to bend over repeatedly, thanks to the two push buttons on the seca 874 that can be turned on and off with a toe tap. Yet another labor-saving feature is the mother/child-function that weighs a child held in its mother's arms. Transport of the scale is made easy by the integrated handle, very low device weight of 8.9 lbs and a carrier bag of feather-light nylon, which is included in the supply schedule.

Intelligent details.

Special mention should be made of the double display that lets the patient and medical personnel read the results from two different directions at the same time. Thanks to a protected toggle switch and a manual power switch, which prevents unintended activation during transport, the function of the scale is always ensured even under difficult conditions.



The flat design and slip-proof surface ensure easy access and a secure stance.

seca 874

Technical data

- Capacity: 440 lbs / 200 kg
- Graduation: 0.1 lbs < 330 lbs > 0.2 lbs / 50 g < 150 kg > 100 g
- Dimensions (WxHxD): 12.6 x 2.4 x 14.2" / 321 x 60 x 361 mm
- Device weight: 8.9 lbs / 4 kg
- Power supply: batteries
- Functions: mother/child-function, auto-HOLD, lbs/kg switch-over, automatic switch-off
- Optional: stadiometer seca 217 combined with adapter element seca 437, carrying case seca 414 (for combination with stadiometer and adapter element) and seca 421



The large leveling bases on the bottom ensure stability on any type of flooring.



The two push buttons on the rear side of the scale can be operated with the foot.

seca 869

Flat scale with cabled remote display



- Separate display and operating element with infinitely adjustable tilt angle for convenient use.
- Large, easy-to-read $\frac{3}{4}$ -inch-high LCD figures.
- BMI function to determine nutritional status.

seca 869:

What was difficult is now easy.

The seca 869 is exactly right for those looking for a scale with practical cabled remote display and a handle for flexible use. This scale has many functions that make life easier. It's simple to use and easy to clean.



Read the display while at your desk.

The seca 869 makes your work easy. Simply place the display on your desk and the scale in front of it. Connect the two with the 7¼ ft. cable and you're ready to work. It's also easy to read the results since the tilt angle of the display is infinitely adjustable.

Stability even on problematic surfaces.

An especially large leveling base gives the seca 869 great stability. The display and operating element is similarly equipped with a stable rubber base and with an easy-to-use self-explanatory keyboard. Storage is also effortless, thanks to its compact dimensions and light weight.

Measuring with all essential functions.

The seca 869 has several convenient functions required in routine medical practice: BMI to determine nutritional status, HOLD to keep the measurement results in the display, TARE to tare the weight, e.g., to weigh a small child held in a parent's arms.



Compact design and low weight for great carrying comfort and universal usage.



Leveling base on the bottom of the scale base for safe and secure placement.



Infinitely adjustable tilt angle in display and operating element for installation on a desk or wall.



seca 869

Technical data

- Capacity: 550 lbs / 250 kg
- Graduation: 0.2 lbs < 330 lbs > 0.5 lbs / 100 g < 150 kg > 200 g
- Dimensions (WxHxD): 12.6 x 2.4 x 14" / 321 x 60 x 356 mm
- Dimensions/platform (WxHxD): 11 x 2.4 x 10.9" / 280 x 60 x 278 mm
- Cable length: 78.7" / 2 m
- Weight: 9.9 lbs / 4.5 kg
- Power supply: Batteries, power adapter optional
- Functions: TARE, HOLD, BMI, lbs/kg switch-over, automatic switch-off
- Optional: stadiometer seca 217 combined with adapter element seca 437, switch-mode power adapter seca 400, carrying case seca 414 (for combination with stadiometer and adapter element) and seca 421