

BetterLiving

PressureSurface DoubleLayer



An excellent multi-use pressure solution for nursing homes, hospitals or home care environments.

This quality two layer pressure surface is ideal for patients with partial mobility. A top layer of high density, temperature sensitive memory foam conforms to the body's shape, working in combination with firm medium density core for maximum support. Foam is covered by an inner cotton cover for added protection.

The polyurethane top cover features multi-stretch material for reduced shear and friction. With waterfall skirt to provide protection from moisture ingress, seams are stitched for strength, then welded or tape sealed for maximum infection control. A PVC slip-resistant base helps maintain product position on bed. The cover is also pre treated for protection against dust mites, mould, mildew and bacteria and includes an asset management panel for recording facility details.

- ✓ **dual layer surface for medium risk**
- ✓ **great multi-purpose solution for general use in any environment**
- ✓ **multi-density foam layers for maximum comfort and support**
- ✓ **max user weight 130 kg**

Formerly known as Deluxe

Specifications detailed over

BetterLiving


PressureSurface DoubleLayer



Vapour permeable, bacteriostatically treated, fire retardant and water resistant cover



Seams are stitched, then welded or tape sealed for maximum infection control

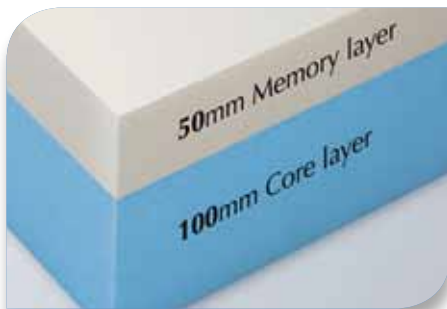

Memory foam combined with medium density core; protected by an inner cotton cover



200 mm extension available for extended sleeping platforms

ITEM	DoubleLayer	Double Layer Extension
CODE	BSPS0452	BSPS0472
CAPACITY	130 kg	NA
HEIGHT	150 mm	150 mm
WIDTH	865 mm	865 mm
LENGTH	1960 mm	200 mm
TOP COVER	Latex free polyurethane	
	6 year body, 2 year cover	
WARRANTY	Product must be removed from vacuum packaging prior to the 'Reconstitute By' date printed on the carton. Failure to do so will void the product warranty.	
TURN	2 way	
COMPLIANCE	ARTG 176837	
FIRE RETARDANCY	Conforms to BS7177:2008 (medium hazard)	
CARE	To clean, wipe down the top cover and PVC base with warm water and a solution of sodium hypochlorite or similar (up to 10,000 ppm available chlorine). Top cover can be machine washed up to 80°C.	
ACCESSORIES	Replacement top cover BSPS0464	
CARTON SIZE	950 x 330 x 330 mm	