

BetterLiving

PressureSurface TripleLayer



The three layered design of this high risk pressure surface is ideal for patients with poor mobility or uncertain skin sensitivity.

Made from premium material, each foam layer is made of varying densities to maximise the redistribution of pressure.

A top layer of high density, temperature sensitive, memory foam conforms to the body's shape, working in combination with the higher resistance support layer to disperse pressure evenly. The special immersion, slow recovery memory foam is ideal for vulnerable areas such as the heels. The core layer and edges provide form, strength and stability. Foam is covered by an inner cotton cover for added protection.

The polyurethane top cover features multi-stretch material for reduced shear and friction. With waterfall skirt to provide protection from moisture ingress, seams are stitched for strength, then welded or tape sealed for maximum infection control. A PVC slip-resistant base helps maintain product position on bed. The cover is also pre treated for protection against dust mites, mould, mildew and bacteria and includes an asset management panel for recording facility details.

- ✓ **six hinge system to conform to profiling beds**
- ✓ **excellent support and pressure redistribution**
- ✓ **three tiered design for people with poor mobility**
- ✓ **max user weight 160 kg**

Formerly known as 3 Layer

Specifications detailed over

BetterLiving


PressureSurface TripleLayer



Vapour permeable, bacteriostatically treated, fire retardant and water resistant cover



Seams are stitched, welded and taped for maximum infection control

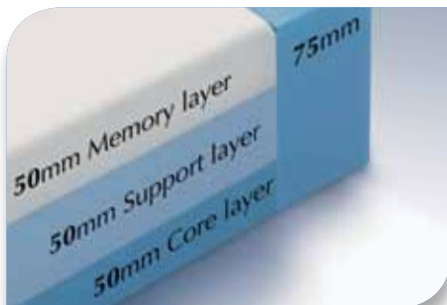

Multi-density foam layers for maximum strength and comfort; protected by an inner cotton cover



Six hinge system to conform to profiling beds

ITEM	TripleLayer	Double Layer Extension
CODE	BSPS0453	BSPS0472
CAPACITY	160 kg	NA
HEIGHT	150 mm	150 mm
WIDTH	865 mm	865 mm
LENGTH	1960 mm	200 mm
TOP COVER	Latex free polyurethane	
	6 year body, 2 year cover	
WARRANTY	Product must be removed from vacuum packaging prior to the 'Reconstitute By' date printed on the carton. Failure to do so will void the product warranty.	
TURN	1 way	
COMPLIANCE	ARTG 176837	
FIRE RETARDANCY	Conforms to BS7177:2008 (medium hazard)	
CARE	To clean, wipe down the top cover and PVC base with warm water and a solution of sodium hypochlorite or similar (up to 10,000 ppm available chlorine). Top cover can be machine washed up to 80°C.	
ACCESSORIES	Replacement top cover BSPS0465	
CARTON SIZE	950 x 330 x 330 mm	