



Etac Aktiv reachers

Award-winning ergonomically designed aids for daily living

 **etac**[®]

Creating Possibilities

Aktiv reachers are light and have ergonomically designed handles to make reaching and gripping easier. The jaws open wide but still grab the smallest item. They are available in a number of models to suit various requirements. All models have a magnet in the handle to help pick up metallic objects.





Aktiv reacher

The jaws of the reacher can be set in two positions. This means that the reacher can be used both horizontally and vertically without having to turn the arm.

The handle is slightly angled and only requires a small hand grip to manoeuvre. This utilises the strength in the hand and provides an effective grip. The entire gripping surface of the jaws is coated with hard-wearing rubber increasing grip-ability and reducing the amount of strength required.



The jaws can be set to two positions.

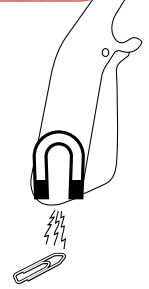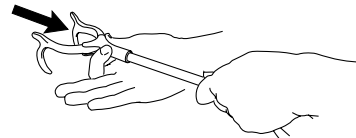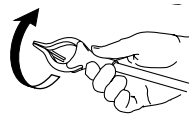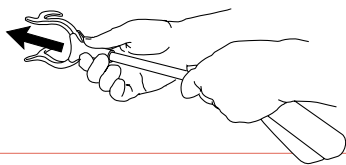


Description	Weight	Item no.	Lifting Capacity
Aktiv reacher, 45 cm 17¾"	111 g 3.9 oz	80505004	The reachers should not be used to lift loads in excess of 2 kg
Aktiv reacher, 70 cm 27½" retail packaging	128 g 4.5 oz	80505103	



Reacher with power grip and hook

This reacher has an 90° angled handle and a forearm support which enables use by people with very limited lifting ability and weak hands. The jaws are provided with hooks to pull or push things away. They can also be set in two positions, which means that the jaws can grip both horizontally or vertically without having to turn the arm.

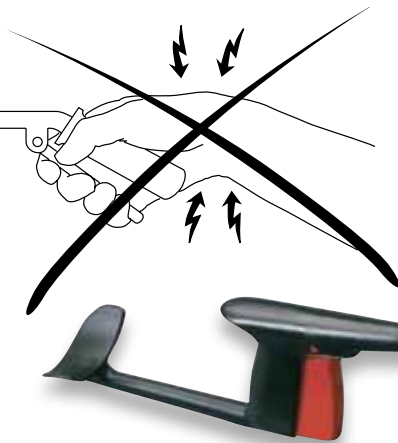

The handle is angled to utilize the strength in the hand and provide an effective grip.



Description	Length	Weight	Order Code	Lifting Capacity
Reacher with Power Grip and Hook	60 cm 23½"	139 g 4,9 oz	80505005	The reachers should not be used to lift loads in excess of 2 kg

Reacher with hook


This handy reacher has many functions. The slightly angled handle utilises the strength in the hand. The jaws are provided with hooks to pull or push things, which also makes it suitable as an aid for dressing and undressing. The jaws can also be rotated and are fitted with patterned rubber on the inside to provide a safe grip.



The hook makes this reacher suitable when dressing and undressing.



Description	Length	Weight	Order Code	Lifting Capacity
Reacher with Hook	45 cm 17¾"	117 g 4.1 oz	80505106	The reachers should not be used to lift loads in excess of 2 kg

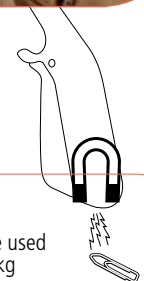



 **etac**[®]
Creating Possibilities

“Everyone should be able to live a free and independent life and pursue their dreams regardless of any physical circumstances.”

 = Standard Stock Item