


# BIOS living

## Ultra Slim Glass Scale

Instruction Manual

SC403

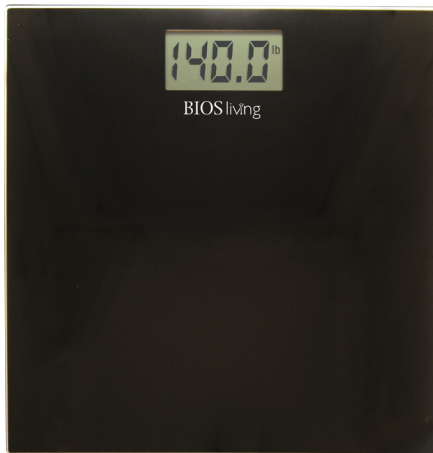

## Introduction

Thank you for purchasing the BIOS Living Ultra Slim Glass Scale. This Scale was developed with state-of-the-art electronics; this precision device has a capacity of 350 lbs / 158 kg, and is equipped with an auto power ON/OFF function, auto zero reset.

## Features

- LCD 3¾" x 1¼" / 9.5 x 3 cm
- Shatterproof glass platform
- Accurate to 350lbs / 158kg
- LB/KG/ST switchable
- Low battery indicator
- Battery included (1 x CR2032)
- 5 year warranty

## Starting your Scale

1. The lithium battery (CR2032) is **included** with your bath scale. Remove the battery cover from the back of the scale, by lifting the black tab. Remove protective film from battery and close the battery cover.
2. Before starting the scale, change the unit of measurement. On the back of the scale there is a small button labeled "KG/LB/ST". Press the button to choose the desired measurement.

## Operation

1. Place the scale on a hard, flat surface (avoid carpet or other soft surfaces). **NOTE:** If the scale is placed on a soft surface, such as a carpet, the results may be inaccurate or an error message may occur.
2. Step onto the scale slowly and gently (do not jump onto the scale). Stand vertically without swaying and place your feet evenly apart.  
**NOTE:** If the display shows "C", this is an indication of an inaccurate reading. If this occurs repeat step 2.
3. Your weight will be displayed on the screen; it will flash three times and then stop at your correct weight.
4. Step off the scale, a new user can use the scale immediately afterwards.
5. The scale will automatically turn off after 8 seconds when not in use.

## Troubleshooting

Error Displayed	Reason for Error
Lo	Indicates low battery
oL	If the scale is overloaded - Over capacity (350 lbs / 158 kg) remove weight immediately from scale
Nothing is displayed on the screen	Replace the battery

## **Advice for Use and Care**

- Avoid using platform when wet - platform will be slippery when wet. Keep scale dry.
- Avoid an extremely high or low temperature environment.
- Avoid excessive impact or vibration to the unit.
- Clean the scale with a damp cloth and avoid water getting inside scale.
- Remove battery when the scale is not being used for long periods of time.
- Store scale in horizontal position.

## **Five Year Warranty**

If this product proves to be defective in material or workmanship within five years of purchase, please return it to the address below. It will be repaired or replaced without charge upon receipt of the unit prepaid with \$5.00 to cover handling, packaging and return postage. Please include proof of purchase, your full name, address, daytime phone number or email address.

Thermor Ltd.  
16975 Leslie Street  
Newmarket, ON L3Y 9A1  
[www.biosmedical.com](http://www.biosmedical.com)