



Leone Bone and Tissue Profiler Tool

Leone Bone profiler tool is used to assist in removing and profiling bone and Soft Tissue before restoration.

Made of Stainless Steel

The Bone and Soft Tissue profiling device aids in uncovering an implant after the healing period and prior to restoration. Three guide pins are available, one for each Leone Implant Connection. The guide pin is installed in the implant attachment socket, and the cutting device slides over the guide pin and is rotated with your implant hand piece, to cut away excess bone and soft tissue that has grown over the top of the implant.


The Cover Screw or Healing cap must be first removed

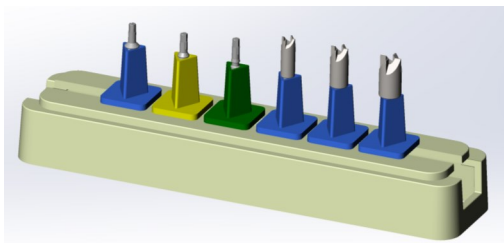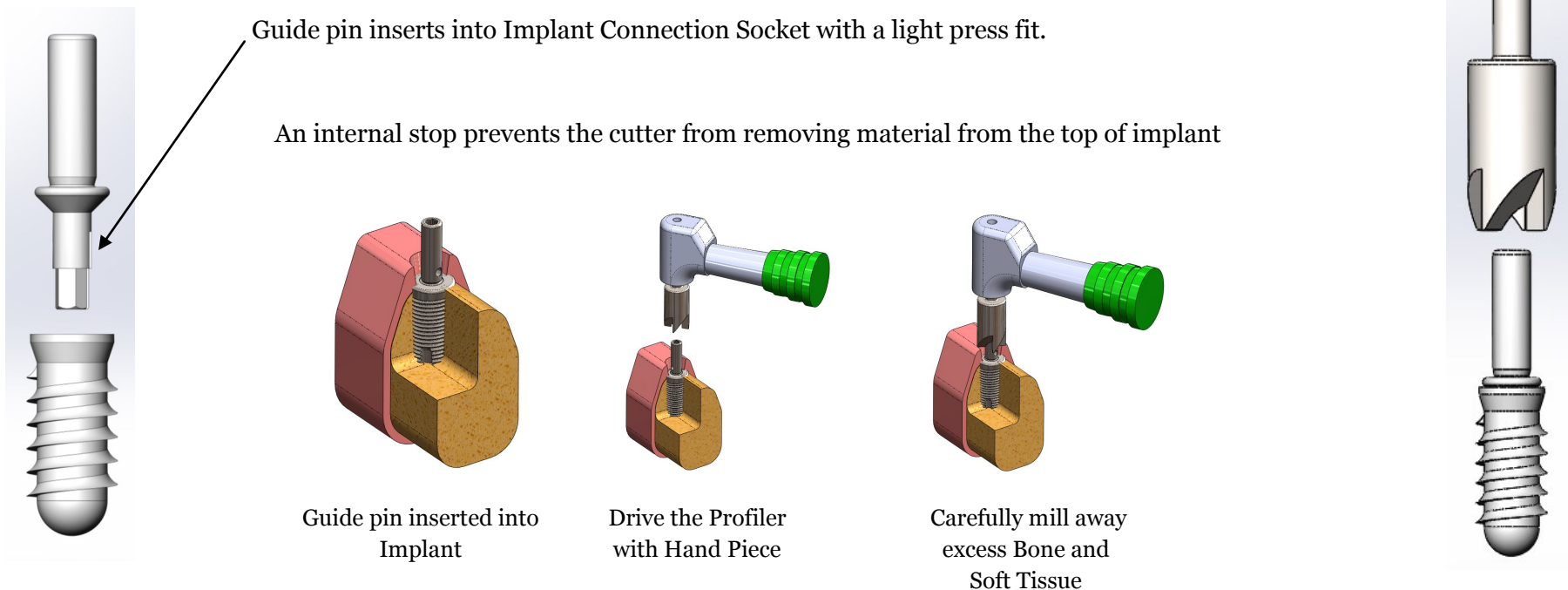

Part Numbers

Guide Pins

- 1156-1500-33
- 1156-1500-41
- 1156-1500-48

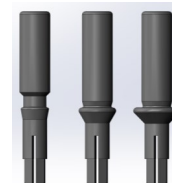

Profilers

- BP100
- BP200
- BP300

Bone Profiler components in kit form... P/N 1156-1500-00