

# the prophy range



The largest and most extensive range of prophy brushes, cups, points and prophy angles on the market today.

premium  
plus™

## Shanks

A quick guide to the advantages of each shank and type of handpiece required.



**CA/RA Metal Shank**  
Most popular style connection



**CA/RA Plastic Shank**  
Versatile and economical shank



**CA/RA Metal Short Shank**  
Designed for greater oral cavity access



**SCREW IN SHANK**  
Provides the most stable vibration free use



**SNAP-ON**  
The absence of a shank delivers the most economical & environmentally friendly cup

Snap-on adaptor not supplied with cups

## Shapes

All our cups are available in regular, soft and extra soft.



### RIBBED & WEBBED CUP

Designed with 4 webs and internal ribs, this cup delivers the optimum amount of prophy paste to help remove stains. The cup also has the ability to "flare" out and therefore reach subgingivally and into interproximal areas.

LATEX AND BPA FREE



### SHORT RIBBED & WEBBED CUP

All the benefits of our standard cups but with its shorter length gives the advantage of greater access in to the oral cavity. When used in combination with a short shank and mini CA handpiece, the total size of the cup and handpiece head is reduced by over 30%. Particularly effective for patients with restricted access in the oral cavity or those sensitive to gag reflex.

LATEX AND BPA FREE



### TURBINE BLADED CUP

Designed with nodules on the exterior surface to optimize interproximal polishing, stain and biofilm removal, and thus requiring less pressure and relieving user fatigue. Turbine blades keep paste in the cup reducing splatter whilst still allowing the cup to flare.

LATEX AND BPA FREE

# Colours

The largest range of colours make identifying the grades of hardness simple, whilst adding a splash of colour to the surgery.

Our prophy cups are available in either single or assorted colour boxes.



White



Light Blue



Pink



Lime



Orange



Lavender



Green



Dark Blue



Yellow




Red

**REGULAR CUP**

**SOFT CUP**


**EXTRA SOFT CUP**

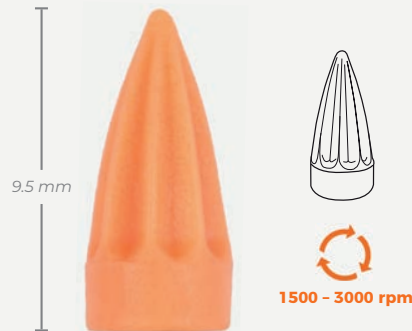


Designed to flare to access interproximal & subgingival areas



Exterior nodules polish interproximal and subgingival areas

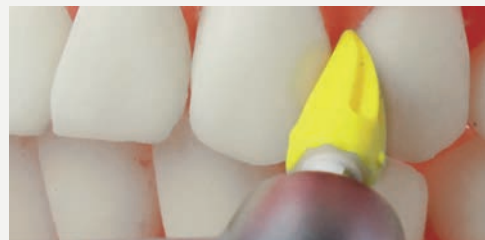

Easy access to tight interproximal areas

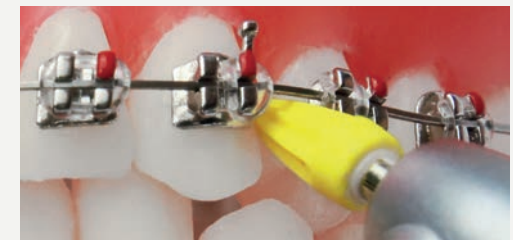

Ideal for orthodontic appliances

**SHORT TURBINE BLADED CUP**

All the benefits of our standard cups but with its shorter length gives the advantage of greater access into the oral cavity. When used in combination with a short shank and mini CA handpiece, the total size of the cup and handpiece head is reduced by over 30%. Particularly effective for patients with restricted access in the oral cavity or those sensitive to gag reflex.

LATEX AND BPA FREE

**POINT - ORTHODONTIC**

Designed for accessing tight or hard to reach areas. Ideal for orthodontic or periodontal procedures where flexibility allows for the most efficient polishing. The fluted ribs also provide ridges which offer the most effective cleaning and polishing surface.

LATEX AND BPA FREE

## Shanks

A quick guide to the advantage of each shank and type of handpiece required.



**CA/RA Metal Shank**  
Most popular style connection



**CA/RA Metal Short Shank**  
Designed to provide greater oral cavity access



**SCREW IN**  
Provides the most stable vibration free use



**SNAP-ON**  
The absence of a shank delivers the most economical & environmentally friendly brush  
Snap on available in flat & tapered only

## Shapes

Which shape or size best suits my needs...

9 mm



1500 - 3000 rpm

### 9MM JUNIOR CUP

Our standard classic shaped nylon bristled brush. The most popular brush for polishing dentition, gold and amalgam fillings.

LATEX AND BPA FREE

7 mm



1500 - 3000 rpm

### 7MM MINI CUP

A smaller version of our popular Junior Cup. All the advantages of the Junior Cup but with greater access in the oral cavity.

LATEX AND BPA FREE



1500 - 3000 rpm

### MICRO BRUSH

Designed for reaching areas with limited access, the nylon bristles are perfect for polishing in and around orthodontic appliances, interdental spaces, implants and cavities.

LATEX AND BPA FREE

# Colours

The largest range of colours make identifying the brush hardness simple and adds a splash of colour to the surgery.

Including Extra Soft which is ideal for work around implants.



White



Light Blue



Clear



Pink



Lime



Orange



Lavender



Green



Dark Blue



Yellow



Red

### REGULAR BRUSH

White available in single colour boxes only.  
Light blue, orange and lavender available in assorted boxes only.

### SOFT BRUSH

Clear available in single colour boxes only.  
Pink, lime, green and dark blue available in assorted boxes only.

### EXTRA SOFT BRUSH



1500 - 3000 rpm

### FLAT BRUSH

Our larger flat brush for greater coverage when polishing dentition, gold and amalgam fillings.

LATEX AND BPA FREE



1500 - 3000 rpm


### TAPERED BRUSH

Our tapered brush gives greater coverage and the ability to pin point delicate areas on implants and orthodontic appliances.


LATEX AND BPA FREE




Junior cup and mini brush are designed for polishing dentition and amalgam fillings

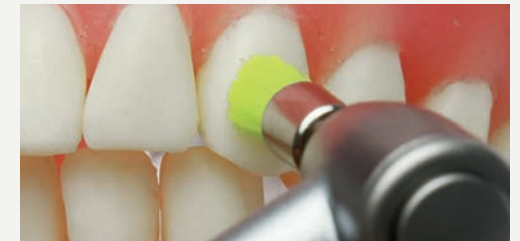

Flat brush easily accesses surfaces and interproximal areas

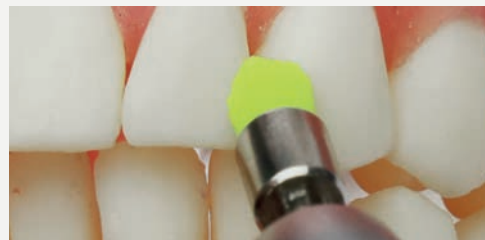

Micro brush designed to access small, hard to reach areas that need precision

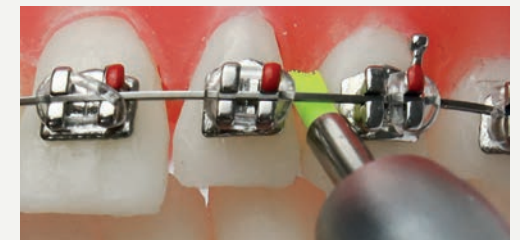

Tapered brush ideal for orthodontic appliances or areas of tight access

## Why use a prophy angle?

Prophy angles fit most standard straight handpieces and this eliminates the need to buy expensive purpose built prophy handpieces. They come individually wrapped to improve infection control and their design prevents the need for disassembling handpieces so as to clean the inner mechanism and gears after use. Prophyl angles save time and this means saving money.

## 360° vs 60° reciprocating

Our standard range of prophy angle cups rotate 360° providing efficient, economical and safe polishing. Our 60° reciprocating heads take efficiency to a new level. They produce no heat and can therefore be used continually, saving time and money. There is no saliva, blood or prophy paste splatter and the lack of motion means no chance of gum or gingiva laceration. Simply a better choice.



## 90° vs 105°

Premium Plus prophy angles are available in two different angled heads: 90° or 105°. Both are ergonomically designed to provide comfortable fatigue free use and have precision crafted gearing to prevent vibration.

The choice is yours.



## 60° Reciprocating

Cups available in singular colours or assorted boxes.

Brushes available in single colour boxes of regular white, soft clear or extra soft yellow. Colours available in assorted boxes only.



Turbine bladed cup holds prophy paste in place but allows the edge to flare producing better contact and more effective polishing. Exterior ridges polish interdentally, saving time and money.

All the advantages of the standard cup with a smaller, shorter cup allowing greater access to all areas in the oral cavity.

Our flat brush for greater coverage when polishing dentition, gold and amalgam fillings.

Our tapered brush gives greater coverage and the ability to pin point delicate areas on implants and orthodontic appliances.

Classically shaped brush provides efficient and effective action for polishing dentition, gold and amalgam fillings.

### REGULAR

### SOFT

### EXTRA SOFT



White Light Blue Orange Lavender Clear Pink Lime Green Dark Blue Yellow Red

# 360° Rotational Clear

Our best selling prophyl angle.

## RIBBED & WEBBED CUP



Designed with four webs and ribs this cup delivers the optimum amount of prophyl paste to help in removing stains.

## SHORT RIBBED & WEBBED CUP



All the advantages of our standard ribbed and webbed cups with a shorter cup allows access to all areas in the oral cavity.

## TURBINE BLADED CUP



Turbo bladed cup holds prophyl paste in place but allows the edge to flare producing better contact and more efficient polishing. Exterior ridges polish interdentially saving time and money.

## REGULAR



Light Blue Lavender

## SOFT



Green Pink

# 360° Rotational Standard

Our most economical prophyl angle. Simple and effective.

## RIBBED & WEBBED CUP



Designed with four webs and ribs this cup delivers the optimum amount of prophyl paste to help in removing stains.

## FLAT BRUSH



Our flat brush for greater coverage when polishing dentition, gold and amalgam fillings.

Brushes available in regular white only.

## TAPERED BRUSH



Our tapered brush gives greater coverage and the ability to pin point delicate areas on implants and orthodontic appliances.

## REGULAR



Lavender White

## SOFT



Green



[www.premiumpluspoland.com](http://www.premiumpluspoland.com)

**European Office**  
**Premium Plus Poland sp. z o.o.**  
ul. Bukowska 27  
62-081  
Wysogotowo  
**Poland**

+48 61 880 1094  
info@premiumpluspl.com  
www.premiumpluspoland.com

**UK Office**  
**Premium Plus UK Ltd.**  
Unit 2, Knighton Heath Ind. Estate,  
Bournemouth,  
Dorset BH11 8NE  
**UK**

+44 1202 611011  
info@premiumplusuk.com  
www.premiumplusuk.com

**Corporate Head Office**  
**Premium Plus International Ltd.**  
Flat 1001, Yuen Long Trading Centre,  
No.33 Wang Yip Street West,  
Yuen Long  
**Hong Kong**

+852 2761 3211  
mailbox@premiumplushk.com  
www.premiumplushk.com

**Manufacturing Facility**  
**Premium Plus China Limited**  
Flat 102, No.122, Tianqiao Road  
Tutang, Changping, Dongguan,  
Guandong  
**China**

+86 769 83756768  
info@premiumpluschina.com  
www.premiumpluschina.com

**USA Office**  
**Premium Plus Dental Supplies Inc.**  
155 Liberty Avenue  
Brooklyn  
New York 11212  
**United States of America**

+1 718 366 8686  
info@premiumplusdentalusa.com  
www.premiumplususa.com

**Japan Office**  
**Premium Plus Japan Co. Ltd.**  
5th Floor, 1-1-31, Nishimidorigaoka,  
Toyonakashi  
Osaka 560-0005  
**Japan**

+81 6 6845 0066  
info@premiumplus.jp  
www.premiumplus.jp