



Creating smiles? You have choices. *Veneer Cementation*

Choice 2 is a light-cured luting cement designed specifically for cementation of porcelain and composite veneers. Choice 2 exhibits superior color stability which is a critical factor in aesthetic veneer cementation. A range of VITA® shades as well as chromatic shades are available allowing for flexibility to achieve high aesthetics. Corresponding try-in pastes (sold separately) are available for shade confirmation prior to final cementation.



- Choice 2 Starter Kit (See reverse for details)
- Choice 2 Try-In Kit
- Refill Syringe (4.5g)
(A1, A2, A3, B1, B3, C2, D2, Translucent, Universal Opaque, Milky Opaque, Milky Bright.)
- Try-In Paste Refill Syringe (2g)


FEATURES

- Light-cured
- Color stable ($\Delta E < 1.2$)
- Ideal film thickness
- Highly filled
- Offered in a multitude of shades
- Corresponding try-in pastes (sold separately)

BENEFITS

- Allows for proper placement of veneers
- Undetectable shade shifting for superior aesthetics
- Ensures veneers are completely seated
- Enhances overall strength and stress resistance
- Veneers will blend well and look natural
- Confirms shade selection prior to cementation


Rockwell Hardness Test

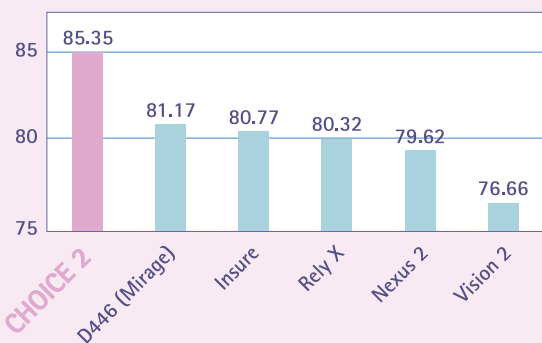

▲ Data from Dr. N. Barghi, University of Texas at San Antonio

Color Stability of Five Luting Resins

$\Delta E < 3$ is not detectable to the human eye

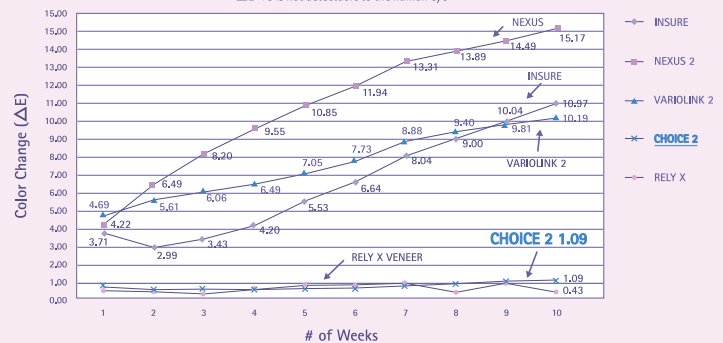



Optimal and Aesthetic Veneers Begin with CHOICE 2

CHOICE 2 STARTER KIT

- 5 Syringes Light-Cured Veneer Cement (4.5g ea.) Shades: 3 x Translucent, 1 x A2, 1 x Milky Bright, 1 Syringe UNI-ETCH® w/ BAC (5g), 1 Syringe 9.5% Porcelain Etchant (5g), 1 Bottle each Bis-Silane™ Parts A & B (6ml ea.), 1 Bottle each ALL-BOND 3™ Parts A & B (6ml ea.), 1 Bottle Porcelain Bonding Resin (6ml), Accessories, Instructions/MSDS.

ENSURES SUCCESSFUL VENEER CEMENTATION



9.5% PORCELAIN ETCHANT 9.5% Buffered Hydrofluoric Acid Gel

- Accurately and efficiently etch a porcelain veneer to ensure optimal bonding.
- Gel material allows accurate placement, full coverage of the veneer, and greater control to ensure the safe usage of the product.



PORCELAIN BONDING RESIN HEMA-Free Unfilled Resin

- A HEMA-free resin designed specifically for bonding porcelain.
- A hydrophobic resin, water will not penetrate the bond thereby eliminating additional shade shifting concerns.



BIS-SILANE™ Two-Part Silane Coupling Agent

- Protects the veneer from contamination and greatly enhances the overall bond of the veneer to the tooth.
- A two-component silane ensures that the material is always fresh, resulting in optimal performance.



ALL-BOND® 3 The Gold Standard Just Got Better

- Easier and faster to use
- More hydrophobic for longer lasting bonds
- Consistent and predictable results



DUO-LINK™ Composite Luting Cement

- Extremely high degree of conversion and stability ensures a long lasting restoration placement
- A universal cement with a proven clinical track record



32% or 37%

BISCO Etchants Polymer Thickened Phosphoric Acid

- Antimicrobial proven effective even after rinsing w/BAC. (Also available w/out BAC)
- Does not leave silica residue thereby optimizing bond strengths
- Superior handling characteristics / stays where it is placed