

	Pediatric	Bicuspid	Bicuspid with deep cervical preps	Bicuspid and small molars	Molars	Molars with deep cervical preps
FXH Series	Composi-Tight® 3D Fusion™ Firm matrix bands with eZ-Place™ technology for easier placement.					
	thickness 0.0015"	NA	FXH100 4.3 mm	FXH150 4.3 mm w/ext 5.9 mm	FXH175 5.3 mm	FXH200 6.1 mm
	height				FXH300 6.1 mm w/ext 7.7 mm	
FX Series	Composi-Tight® 3D Fusion™ Full Curve Bands have non-stick coating and an occlusal placement tab.					
	thickness 0.0016"	NA	FX100 4.0 mm w/ext 4.4 mm	FX150 4.1 mm w/ext 6.0 mm	FX175 5.0 mm w/ext 5.6 mm	FX200 5.4 mm w/ext 6.6 mm
	height				FX300 6.0 mm w/ext 8.7 mm	
SXR Series	Slick Bands™ XR Matrixes have non-stick coating and Grab-Tabs™ for easier placement.					
	thickness 0.0016"	NA	SXR100 4.6 mm	SXR150 3.8 mm w/ext 5.4 mm	SXR175 5.4 mm	SXR200 6.4 mm
	height				SXR300 6.4 mm w/ext 8.6 mm	
SM Series	Slick Bands™ Matrixes have non-stick coating.					
	thickness 0.0016"	NA	SM050 3.2 mm w/ext 4.2 mm	SM100 4.6 mm	SM150 3.8 mm w/ext 5.4 mm	SM175 5.5 mm
	height				SM200 6.4 mm	SM300 6.4 mm w/ext 8.6 mm
M Series	Composi-Tight® bands are contoured in 3 dimensions to provide an accurate contact at the tooth's natural height of contour.					
	thickness 0.0015"	NA	M050 3.2 mm w/ext 4.2 mm	M100 4.6 mm	M150 3.8 mm w/ext 5.4 mm	M175 5.5 mm
	height				M200 6.4 mm	M300 6.4 mm w/ext 8.6 mm
AU Series	The Gold bands are the same as the Composi-Tight® bands.					
	thickness 0.0015"	NA	AU050 3.2 mm w/ext 4.2 mm	AU100 4.6 mm	AU150 3.8 mm w/ext 5.4 mm	NA
	height				AU200 6.4 mm	AU300 6.4 mm w/ext 8.6 mm
B Series	The original bands have a slightly shorter length.					
	thickness 0.0015"	NA	B050 3.2 mm w/ext 4.2 mm	B100 4.6 mm	B150 3.8 mm w/ext 5.4 mm	NA
	height				B200 6.4 mm	B300 6.4 mm w/ext 8.6 mm
MC Series	Clear bands are transparent and translucent to allow for trans-enamel polymerization without interference from metal bands.					
	thickness 0.0020"	NA	MC100 5.0 mm	NA	MC175 5.9 mm	MC200 6.8 mm
	height					

Item
FXH100
FXH150
FXH175
FXH200
FXH300

Item
FX100
FX150
FX175
FX200
FX300

Item
SXR100
SXR150
SXR175
SXR200
SXR300


Item
SM050
SM100
SM150
SM175
SM200
SM300

Item
M050
M100
M150
M175
M200
M300

Item
AU050
AU100
AU150
AU200
AU300

Item
B050
B100
B150
B200
B300

Item
MC100
MC175
MC200

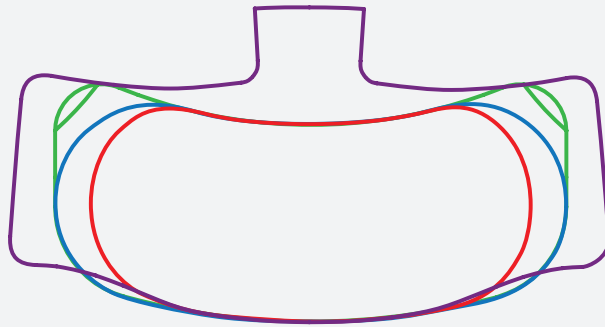

FX Series
SXR Series
M – SM – AU Series
B Series

FX Series
(coloured/nonstick band)

- The longest band (buccal to lingual). Ideal for a large preps, or when using the wide green ring.
- This curved band fits the anatomy of the tooth, making it easier to go subgingivally.
- Includes a grab tab for easy placement.

FX, SXR, SM series
(coloured/nonstick band)

- These bands all have a non-stick coating, making them much easier to remove! (The chemistry in newer dental adhesives and composites can make it challenging to remove non-coated bands.)
- Colour coding makes it easier for dental team members to identify the bands clinically and when reordering.

FXH series
(silver - not pictured above)

- The newest addition! An extra firm matrix band.
- It is the same thickness as the FX but constructed from a denser material, giving it the feeling of being firm.
- Ideal for tight interproximal spaces where other bands may crumple more easily upon insertion.
- Silver colour, curved, and slightly shorter (buccal to lingual) than the FX series.
- Includes a grab tab for easy placement.

B series
(silver)

- The original sectional matrix bands!
- Shortest in length - ideal for Pedo

M (silver), **AU** (gold) and
SM (coloured/nonstick)
series

- Band size & thickness is the same for all.
- Available in Silver (M), Gold (AU), Coloured/Non-stick (SM)